Laparoscopic Instruments

Product Catalog



Aesculap Surgical Technologies - Endoscopy

AESCULAP

Table of Contents

INDEX

e villani	INTRODUCTION	4
	History	4
	Surgical Quality Management	7
AND		
500	INSTRUMENTS	11
131	Scissors, Graspers, Dissectors	12-27
25	Modular System	12-23
-8	Reposable Scissors	27
	Monopolar Electrodes	28
	Suturing	29
	Bipolar Instruments	30-31
	Specialty-Locking Graspers	32
	Specialty Instrument Sets - Advanced Procedure	33-36
	Specially instrument sets - Advanced Procedure	33-30
	OLID ADDITEDC	07.40
74 / /	CLIP APPLIERS	37-40
-	SLS Single-Fire Ligation System	38
-	Challenger Ti-P Multi-Fire Ligation System	39
-	Vascular Clip Applier/Remover	40
1	ENDOSCOPES AND ACCESSORIES	41-43
2 3 7 7 7	HD Endoscopes	42
111011	Unitrac® Retraction & Holding System	43-44
	TROCARS	45-51
The state of the s	Reusable Trocar Systems	46-50
200	Flexible Trocars	51
•		
	VIDEO SYSTEMS	53-55
Walter Com	Light Sources	54
	Light Cables	55
	MISCELLANEOUS	57-62
The second		
A STATE OF THE PARTY OF THE PAR		

History

Blending Old World Craftsmanship with Modern Technology



Aesculap.

Aesculap's success derives from

As the world's largest surgical instrument manufacturer, Aesculap designs and produces surgical instruments, implants and motors for several surgical specialties including: ENT, plastic and reconstructive, thoracic, cardiovascular, orthopedic, gynecologic and general surgery as well as neurosurgery and minimally-invasive surgery.

The company was founded in 1867 by master craftsman Gottfried Jetter in Tuttlingen, Germany. For more than 140 years, the name Aesculap has been synonymous with quality throughout the world.

Today, modern technology blends with Aesculap Old World tradition to provide healthcare professionals with the widest array of surgical instrumentation available from a single manufacturer. Quality craftsmanship

a history of:

- Surgeons' designs united with instrument expertise
- Total commitment to customer service
- Intense communication with our customers
- National sales force
- National hands-on workshops, specialized by discipline
- Decades of international training at the Aesculap Academy – Tuttlingen, Germany

Aesculap's product range:

- Endoscopy
- Surgical instruments for all disciplines

Gottfried Jetter

- SterilContainers[™]
- Neuroendoscopy
- Aneurysm clips
- Power systems
- Ultrasonic aspiration
- Spinal implants
- Orthopaedics
- Per-procedure programs
- Repair services
- Consulting services

U.S. Headquarters



U.S. Headquarters

Our U.S. Headquarters, located in Center Valley, Pennsylvania, was established in 1977 to serve the needs of the U.S. market.

Aesculap's Number One Core Value is simple: *CUSTOMERS COME FIRST.*

Our primary goal is to exceed the expectations of the medical professionals who depend on our products to deliver healthcare to their patients.

Our Customer Service Professionals are dedicated to serving the needs of all of our valued customers.

Please call customer service for:

- Placing orders
- Product information
- Pricing
- Custom product information

Aesculap, Inc. 3773 Corporate Parkway Center Valley, PA 18034 Phone 800-258-1946 Fax 610-791-6880 www.aesculapusa.com

Customer Service:

Phone: 800-282-9000 *Fax:* 610-791-6886

The Worldwide Learning Experience







Aesculap Academy

The Aesculap Academy was established in 1995 to offer surgeons and medical experts the opportunity to exchange ideas about medical issues and discuss the latest techniques in theory and practice.

Designed as a Forum for Communication in Medicine, the Aesculap Academy has dedicated itself to facilitate an interdisciplinary dialogue to provide knowledge transfer from person to person.

Topics range from state-of-the-art medical practices to public health and scientific controversial issues.

Partnerships and associations with internationally recognized institutions guarantee a global focus on current medical topics.

The Aesculap Academy exists in more than 29 countries worldwide and will continue to grow. The Academy headquarters, or Aesculapium, is located in Tuttlingen Germany, and its U.S. facility is located in Center Valley, Pennsylvania.

For more information:

Call:

1-800-258-1946, Ext. 5223

or visit our website: www.aesculap-academy.com



Our Solutions - Your Success





Think Solutions. Think Aesculap.

Through Intelligent Quality Solutions we are able to share the expertise we've gained through decades of experience. It is our tradition of excellence and our commitment to sharing knowledge that has made Aesculap the partner of choice for top hospitals throughout the world. Let Aesculap help design an individual solution to meet your instrument management needs in the SPD and the OR.

Think Solutions. Think Aesculap.

Consulting Solutions

Technical Solutions

IT Solutions

Supply Solutions

Consulting Solutions

Quick Scan Analysis

An onsite assessment of your OR and Sterile Processing Department (SPD) to customize a solution that meets your strategic goals

- Offering you strategic recommendations for long-term process improvements
- Hands-on best practices and benchmark data provide immediate impact and real results

Set Optimization Consulting

Having the right instruments in the right set for the right procedure

- Having only the instruments you most commonly use in a set to minimize reprocessing costs
- Reduced weight of instrument trays and better ergonomics result in less work for the SPD staff
- Inventory can be reallocated to reduce your need to purchase additional back-up supplies





For more information:

Contact your Aesculap Sales Representative for more information on how these asset management programs can be implemented in your facility.

Inventory Analysis

A detailed analysis of your existing surgical instrument inventory by our team of technicians, nurses and instrument experts

- Analyze and document the quality of your current instrument inventory
- Supply recommendations for care and handling that will prolong your instrument life
- Analyze instrument inventory levels to determine if sufficient to meet hospital procedure volume

Technical Solutions



Aesculap Technical Service (ATS)

is committed to ensuring quality performance throughout the life of every instrument. Aesculap's commitment to quality does not end after an instrument is purchased. With instrument repair and refurbishment services, you can be sure that every instrument functions properly and is available when you need it.

The highest performance is maintained through regular service and repair, keeping surgeons satisfied and improving patient care.

Skilled craftsman and trained instrument makers will refurbish your surgical instruments to "like new" condition, quickly and conveniently. Aesculap uses the same precision equipment in our repair facility that is used to manufacture these instruments in Tuttlingen, Germany. Factory spare parts give Aesculap the ability to repair your instrumentation with quality materials that match the original. This protects your instrument from contamination and corrosion, further reducing repair, replacement costs and down time.

ATS can repair general and endoscopic instruments, rigid and flexible endoscopes, and motors. For your convenience, we also service instruments from all major manufacturers.

For repairs and service, please contact:

Aesculap Technical Service - ATS

St. Louis Distribution & Technical Service Center 615 Lambert Point Drive Hazelwood, MO 63042-2609 800-214-3392



Full Service Repairs

- General instruments
- Endo instruments
- Powered instruments
- Endoscopes
- Micro-Neuro instruments
- Cardiovascular instruments
- Visual systems
- SterilContainer® system

Enhancements

- Gold dipping
- Diamond dusting
- Color coating
- Custom modifications
- Etching
- Re-insulation

Refurbishment

- Single instruments
- Complete sets

Capitated and Risk Sharing Programs

- Reduce operating expenses
- Improve resource management
- Simplify the budgeting process





For more information, please contact ATS at 800-214-3392

IT Solutions



Instrument Management System

Instacount[™] instrument tracking system puts critical information at your fingertips.

instacount PLUS is a comprehensive instrument management system that addresses the asset management needs of sterile processing departments. The use of barcode technology allows for easy, accurate and immediate execution of many workflow processes. The software is designed to reside on a hospital server and interface with a network of computers. Comprehensive reporting options allow managers to make critical business decisions based on accurate, real-time data.


- Tracks location of instrument set/item
- Verifies sterilization type/parameters and load tracking
- Monitors staff productivity

UNICOS/Data Matrix

Tracking your supplies down to the instrument level

- Traceability of instruments down to the individual patient
- Machine readable technology saves your technicians time
- Tracking of specialty one-of-a-kind or expensive instruments reduce loss
- All Aesculap instruments, size permitting, come with a data matrix label at no additional cost



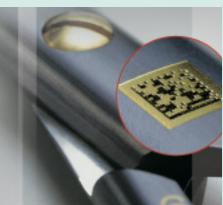





For more information:

Please contact your Aesculap representative for further details or call Customer Service at 800–282–9000

Supply Solutions



PCMS

Aesculap provides high quality surgical products and services to a facility on a fee per procedure basis.

- Provide new technology to your surgeons quickly and without capital investment
- Not only do you receive new instruments, but also the services that go towards maintaining the instruments for the life of your contract
- Optional on-site technician will process and assemble the instrument sets
- Ensures sets are ready every time, removing the worry of instrument readiness

Our sales representatives work with each surgical and processing team to analyze existing sets, procedure volumes, and costs. *PCMS* guarantees optimal results with high quality, professionally maintained instruments.





ReadyCaseSM

The ReadyUseSM is an innovative surgical supply and instrument delivery system that will change the way your OR operates. This comprehensive service includes daily delivery of instrument sets, reusable surgical packs, disposable accessory packs and single sterile items on a per-procedure basis. The sets are customized to meet individual surgical needs.

ReadyCaseSM can help your facility solve hospital Operating Room and Central Sterile issues:

- Significantly reduce procedure costs
- Eliminate repair and maintenance of instruments
- Decrease in-house sterile processing demands
- Reduce operating room set up and turnaround time
- Obtain quality surgical instrumentation without capital expense

For more information:

Please contact your Aesculap representative for further details or call Customer Service at 800–282–9000

ReadyCaseSM is a service mark of SRI Surgical.





QUICK-SNAP™ Laparoscopic Instruments


Aesculap's SOVEREIGN® line of instruments for minimally invasive surgery combines precision craftsmanship with modern materials and design.

Optional Handles for Customizing Instruments - Interchangeable for 5 mm and 10 mm instruments







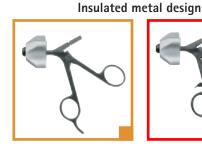





PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet)

PM953R (In-Line)

			COMPONE	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	P0001R				
	Hook Scissors, 5 mm, 31 cm	P0600R	PM973R		
	DOOGED				
	PO025R Hook Scissors, 5 mm, 42 cm	P0744R	PM981R		
	HOOK SCISSOIS, 5 IIIII, 42 CIII	FU/44N	FINISOLU		
1/1	P0003R				
	Micro-Dissecting Scissors, 5 mm, 31 cm	P0602R	PM973R		
	Order by component				
	Micro-Dissecting Scissors, 5 mm, 42 cm	P0745R	PM981R		
	· · · · · · · · · · · · · · · · · · ·			_	
1/1	PO004R				
	Metzenbaum Scissors, serrated,	P0603R	PM973R		
6	Tungsten Carbide inserts, 5 mm, 31 cm	PU603K	PIVI9/3K	_	_
	P0027R				
	Metzenbaum Scissors, serrated,				_
	Tungsten Carbide inserts, 5 mm, 42 cm	P0746R	PM981R		
1/1	PM015R				
	Metzenbaum Scissors, insulated to tip,				
	Tungsten Carbide Inserts, 5 mm, 31 cm	PM697R	PM973R		

SOVEREIGN® QUICK-SNAP™ design











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)


PM951R (Ratchet)

PM953R (In-Line)

			COMPONI	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO090R Premium Metzenbaum Scissors, 5 mm, 31 cm	P0790R	PM973R		
1/1	PO002R Peritoneal Scissors, serrated, 5 mm, 31 cm	P0601R	PM973R	•	
1/1	PO006R Peritoneal Scissors, serrated, 5 mm, 31 cm	P0625R	PM973R		
1/1	PO011R Baby Metzenbaum Scissors, serrated, Tungsten Carbide inserts, 5 mm, 31 cm PO028R	P0655R	PM973R	٠	
	Baby Metzenbaum Scissors, serrated, Tungsten Carbide inserts, 5 mm, 42 cm	P0754R	PM981R		•
1/1	PM016R Baby Metzenbaum Scissors, insulated to tip, Tungsten Carbide inserts, 5 mm, 31 cm	PM698R	PM973R		
1/1	PO091R Premium Baby Metzenbaum Scissors, 5 mm, 31 cm	P0791R	PM973R		





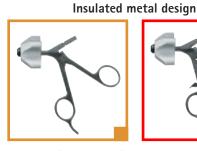





PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet) PM953R (In-Line)

			COMPONENTS		
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO022R CROCE/OLMI Dissecting and Suture Scissors with Knot Pusher, 5 mm, 31 cm	P0727R	PM973R	÷	÷
1/1	PO550R Durogrip TC Needle Holder, straight, 5mm, 31cm	PO628R	PM973R	÷	
1/1	PO559R Durogrip TC Needle Holder, curved, 5mm, 31cm	P0668R	PM973R		**
1/1	PO101R Atraumatic Fixation Forceps, serrated, 5 mm, 31 cm	P0607R	PM973R		•
annua di sa	PO199R Atraumatic Fixation Forceps, serrated, 5 mm, 42 cm	P0747R	PM981R	•	•
1/1	PO129R Tubal Fixation Forceps, 5 mm, 31 cm	P0645R	PM973R		
	PO200R Tubal Fixation Forceps, 5 mm, 42 cm	P0748R	PM981R		
1/1	PO210R Oval Grasper, 5 mm, 31 cm	P0763R	PM973R		
100	PO207R Oval Grasper, 5 mm, 42 cm	P0760R	PM981R		

SOVEREIGN® QUICK-SNAP™ design











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet) PM951R (Ratchet)


Ratchet) PM953R (In-Line)

			COMPON	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO211R				
	Curved Grasper, 5 mm, 31 cm	P0764R	PM973R		
CHRISTIAL	PO208R			_	
369	Curved Grasper, 5 mm, 42 cm	P0761R	PM981R		
1/1	PO149R	DOCOOD	DMOZOD		
	Fixation Forceps, fenestrated, 5 mm, 31 cm	P0680R	PM973R		
	PO206R				
	Fixation Forceps, fenestrated, serrated, 5 mm, 31 cm	P0758R	PM973R		
1/1					
	PO100R				
	Atraumatic Dissecting Forceps, fine tip, 5 mm, 31 cm	P0605R	PM973R		
and the same of th	3 11111, 31 C111	1 000311	1 10137 311	_	
1/1	PO102R				
	Maryland Dissecting Forceps,	P0608R	PM973R		
annull b	5 mm, 31 cm	PU608K	PIVI973K		
	PO202R Maryland Dissecting Forceps, 5 mm, 42 cm	P0752R	PM981R		
1/1	Maryianu Dissecting Forceps, 5 mm, 42 cm	FU/52N	FINISOLU	_	
1/1	PM165R				
The state of the s	Maryland Dissecting Forceps,	D14D		_	
	insulated to tip, 5 mm, 31 cm	PM699R	PM973R	_	
1/1	P0218R				
	Crile Dissecting Forceps, Curved,			_	
	5 mm, 31 cm	P0771R	PM973R		
1	PO213R				
	Crile Dissecting Forceps, Curved,	P0766R	PM976R		
	5 mm, 37 cm	FU/66K	FIVI3/6K		



PO958R (Non-Ratchet)



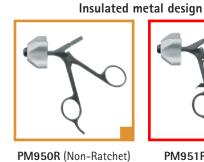






			COMPONE	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO141R Overholt 60° Dissecting Forceps, 10 mm, 31 cm	P0661R	PM975R		::
1/1	PO103R Right Angle, 90° Dissecting Forceps, 10 mm, 31 cm	P0606R	PM975R		**
1/1	PO231R Mixter Dissector, Right Angle, 90° 5 mm, 31 cm	P0784R	PM973R	÷	•
1/1	PO187R Right Angle, 90° Dissecting Forceps, 5 mm, 31 cm	P0730R	PM973R		**
1/1	PO147R Dolphin Nose Dissecting Forceps, 5 mm, 31 cm	P0678R	PM973R	•	::
1/1	PO151R Duck Bill Dissecting Forceps, 5 mm, 31 cm	P0667R	PM973R		::

SOVEREIGN® QUICK-SNAP™ design





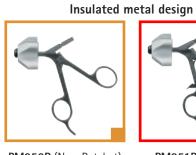





PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet)

PM953R (In-Line)

			COMPONI	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO214R Duck Bill Blunt Grasper, 5 mm, 37 cm	P0767R	PM976R		
1/1	PO190R Bullet Nose Grasper, 5 mm, 31 cm PO203R	P0735R	PM973R	•	
	Bullet Nose Grasper, 5 mm, 42 cm	P0753R	PM981R		
1/1	PO107R Aggressive Grasping Forceps, 5 mm, 31 cm	P0613R	PM973R		•
1/1	PO109R Alligator Grasping Forceps, 5 mm, 31 cm	P0615R	PM973R		•
1/1	PO216R Bowel Grasping Forceps, 20 mm Jaw, 5 mm, 31 cm	P0769R	PM973R	•	•
- 100	PO217R Bowel Grasping Forceps, 20 mm Jaw, 5 mm, 42 cm	P0770R	PM981R	•	**
1/1	PO112R Dorsey Intestinal Forceps, 40 mm Jaw, 5 mm, 31 cm PO204R	P0618R	PM973R		**
The second secon	Dorsey Intestinal Forceps, 40 mm Jaw, 5 mm, 42 cm	P0755R	PM981R		











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet) PM953R (In-Line)

		_	COMPONENTS		
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO191R Universal Grasper, Fenestrated 5 mm, 31 cm	P0736R	PM973R	÷	**
1/2	PO215R Glassman Clamping Forceps, 60 mm Jaw, 10 mm, 31 cm	P0768R	PM974R	÷	**
1/2	PO118R DeBakey Atraumatic Forceps, 50 mm Jaw, 10 mm, 37 cm	P0624R	PM977R	÷	::
1/1	PO137R DeBakey Atraumatic Forceps, 50 mm Jaw, 5 mm, 42 cm	P0657R	PM981R	÷	::
1/1	PO113R Babcock Grasping Forceps, Single Action, 5 mm, 31 cm	P0619R	PM973R	÷	::
1/1	PO212R Babcock Grasping Forceps Double Action, 5 mm, 31 cm PO209R	P0765R	PM973R		
	Babcock Grasping Forceps, Double Action, 5 mm, 42 cm	P0762R	PM981R		

SOVEREIGN® QUICK-SNAP™ design











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet)

PM953R (In-Line)

			COMPONE	ENTS	
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO114R Babcock Grasping Forceps, Single Action, 10 mm, 31 cm PO205R	P0620R	PM974R		
	Babcock Grasping Forceps, Single Action, 10 mm, 42 cm	P0757R	PM982R		
1/2	PO196R Babcock Grasping Forceps, Double Action, 10 mm, 37 cm	P0741R	PM977R		::
1/1	PO117R Claw Grasper, 2 x 3 teeth, 10 mm, 31 cm	P0623R	PM974R		**
1/1	PO104R Allis Atraumatic Grasping Forceps, 5 mm, 31 cm	P0610R	PM973R		::
1/1	PO111R Allis Grasping Forceps, 5 mm, 31 cm	PO617R	PM973R		**
1/1	Universal Grasping Forceps, Fenestrated 5 mm, 31 cm	P0647R	PM973R		•
19111111	PO229R Universal Grasping Forceps, Fenestrated 5 mm, 42 cm	P0783R	PM981R	\mathbf{r}	•







Insulated metal design

PM950R (Non-Ratchet)

PM951F





PO958R (Non-Ratchet)

PO959R (Ratchet)

PM951R (Ratchet)

PM953R (In-Line)

		COMPONENTS			
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES
1/1	PO132R Grasping Forceps, 2 x 4 teeth, serrated, 5 mm, 31 cm	P0648R	PM973R		-
Evenous -	PO201R Grasping Forceps, 2 x 4 teeth, serrated, 5 mm, 42 cm	P0749R	PM981R		•
1/1	PO134R Grasping Forceps, fine-toothed, non-serrated, 5 mm, 31 cm	P0650R	PM973R		•
1/1	PO105R Heywood-Smith Tissue Forceps, fenestrated, 10 mm, 31 cm	P0611R	PM975R	÷	::
1/1	PO130R Atraumatic Grasping Forceps, fenestrated, 5 mm, 31 cm	P0646R	PM973R		::
1/1	PO221R Fundus Grasper, 5 mm, 31 cm, PO230R	P0774R	PM973R		-
	Fundus Grasper, 5 mm, 42 cm,	P0782R	PM981R		
1/1	PO222R Short Fenestrated Grasper, 5 mm, 31 cm PO223R	P0775R	PM973R		
	Short Fenestrated Grasper, 5 mm, 42 cm	P0776R	PM981R		

SOVEREIGN® QUICK-SNAP™ design











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet)

PM953R (In-Line)

			COMPONENTS			
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES	
1/1	PO224R Tenaculum Forceps 5 mm, 31 cm	P0777R	PM973R			
	PO225R Tenaculum Forceps 10 mm, 31 cm	P0778R	PM974R			
1/1	PO500R Large Biopsy Spoon Forceps, 10 mm, 31 cm	P0705R	PM974R		**	
1/1	PO501R Biopsy Forceps, Cup, fenestrated, 5 mm, 31 cm	P0626R	PM973R			
	PO509R Biopsy Forceps, Cup, fenestrated, 5 mm, 42 cm	P0750R	PM981R			
1/1	PO502R Biopsy Forceps, Box, 5 mm, 31 cm	P0627R	PM973R			
	PO512R Biopsy Forceps, Box, 5 mm, 42 cm	P0759R	PM981R			
1/1	PO505R Biopsy Forceps, with teeth, 5 mm, 31 cm	P0652R	PM973R			
	PO511R Biopsy Forceps, with teeth, 5 mm, 42 cm	P0756R	PM981R			
1/1	PO510R Biopsy Forceps, Hook, 5 mm, 42 cm	P0751R	PM981R			











PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)


PM951R (Ratchet) PM953R (In-Line)

	ORDER BY COMPONENT	JAW INSERT	COMPONENTS INSULATED OUTER TUBES	OPTIONAL HANDLES
1/1	Atraumatic Parenchyma Fixation Forceps, 5 mm, 31 cm	PM609R	PM973R	••
1/1	Dissecting Forceps, delicate, 5 mm, 31 cm	PM614R	PM973R	••
1/1	Pike Mouth Aggressive Grasper, 10 mm, 31 cm	PM621R	PM974R	••
1/1	Large Metzenbaum Scissors, TC, serrated 10 mm, 31 cm	PM604R	PM974R	•

SOVEREIGN® QUICK-SNAP™ design





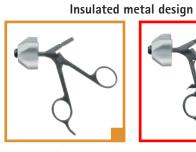





PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)


PM951R (Ratchet)

PM953R (In-Line)


			COMPONENTS		
	ORDER BY COMPONENT	JAW INSERT	INSULATED OUTER TUBES	OPTIONAL HANDLES	
1/2	Endo Lymph Node Fixation Forceps, 10 mm, 31 cm	PM688R	PM974R	••	
1/2	Pennington Tissue Grasping Forceps, 10 mm, 37 cm	PM679R	PM977R	••	
1/2	Duval Atraumatic Fixation Forceps, 10 mm, 37 cm	PM612R	PM977R	••	
1/2	DeBakey Atraumatic Forceps, 40 mm Jaw, 5 mm, 31 cm	PM657R	PM973R	••	

Aesculap Sovereign® Laparoscopic Instruments Instructions for quick and easy assembly and disassembly For use with PM950R and PM951R, insulated metal handles

ASSEMBLY

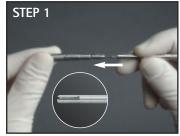


Slide the inner metal tube onto the jaw insert with its four spring tongues towards the working end until the tongues snap into place.

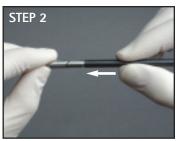

Slide the outer insulation tube as far as it will go over the inner metal tube.



Hold the shaft firmly by its tip and pull back the handle rotation knob with fingers and maintain hold. Fully insert the shaft and release rotation knob.

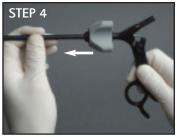


Test that the instrument is secured in the handle by pulling on the shaft. If shaft pulls out of handle, then repeat step 3 and step 4.



Connect handle with jaw insert ball joint by depressing the black buttons on both sides of handle, continue to depress and lift rear ring of handle under ball joint of jaw insert until captured by handle. Release buttons.




Test instrument by opening and closing handle to verify that the jaw has been inserted properly.

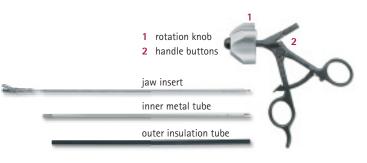


DISASSEMBLY



Depress black buttons on both sides of handle, continue to depress, pull rear ring down to release ball of jaw insert. Release buttons.

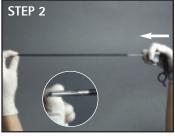

Pull back rotation knob and pull out shaft with other hand. Release rotation knob

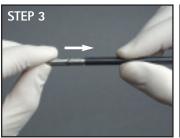


Hold the jaw insert at the tip and remove the outer insulation tube.




Hold the jaw insert at the tip and remove the inner metal tube.

Aesculap Sovereign® Laparoscopic Instruments Instructions for quick and easy assembly and disassembly For use with QUICK-SNAP™ PO958R and PO959R handles


ASSEMBLY

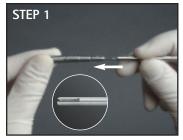


Slide the inner metal tube onto the jaw insert with its four spring tongues towards the working end until the tongues snap into place.

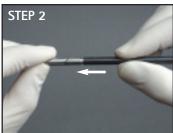

Slide the outer insulation tube as far as it will go over the inner metal tube.

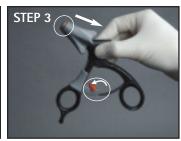

If the red ring is not visible, pull the rotation knob back. If you are using the ratchet style handle, deactivate ratchet by flipping down the red toggle switch.

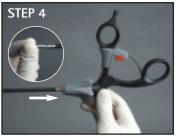


Hold the shaft firmly by its tip and the handle behind the rotation knob. Push the shaft into the handle until it locks in place and the red ring disappears.

Note: Do not hold onto the movable handle component.




Test the assembly of the instrument by opening and closing the handle to ensure the jaw is opening and closing.

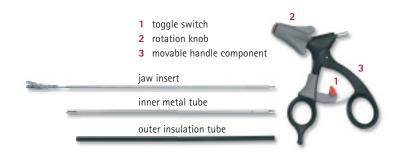


DISASSEMBLY



If ratchet style handle, deactivate ratchet by flipping down red toggle switch.

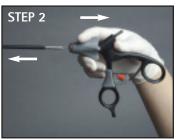

Pull rotation knob back until it locks into place. A red ring will become visible. Once visible remove the shaft by pulling the instrument shaft out of the handle.

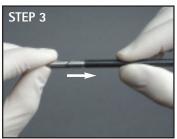

Hold the jaw insert at the tip and remove the outer insulation tube.



Hold the jaw insert at the tip and remove the inner metal tube.

MODULAR INSTRUMENT SPARE PARTS

Insulated Outer Tubes



PM973R 5 mm, 31 cm PM976R 5 mm, 37 cm PM981R 5 mm, 42 cm



PM974R 10 mm, 31 cm **PM982R** 10 mm, 42 cm



PM975R 5/10 mm, 31 cm PM978R 5/10 mm, 37 cm



PM977R 10 mm, 37 cm

P0950400	Cautery Pin for PO handles, 5 pkg.
PL170251	Cautery Pin for PM handles, 5 pkg.
PM950250	Rotation Knobs for PM instruments, grey, 5 pkg.
GK246	Disposable Monopolar Cable, 12 ft.
GK245	Reusable Monopolar Cable, ERBE T 4 mm diam., 12 ft.
US354	Reusable Monopolar Cable, VALLEYLAB, CONMED, BIRCHER, BOVIE, ERBE ICC, 4 mm diam., 10 ft.
GK469R	Cleaning Brush, 3.7 mm diam. for Inner Tube and Outer Tube, 5 mm diam., 12 pkg.
PM995R	Cleaning Brush, 10 mm diam., 6 pkg.
PM600201	Inner Metal Tube, 31 cm
PM618201	Inner Metal Tube, 37 cm
PM434201	Inner Metal Tube, 42 cm

Reposable Scissors

SOVEREIGN® Reposable Scissors



• Reusable Handle and Disposable Jaw Insert

- Sharp scissors tip EVERY TIME with disposable jaw insert
- Cost effective alternative compared to 100% disposable scissors
- Three different scissors options
- Convenient handle choices

•		
	ITEM NO.	DESCRIPTION
1/1	P0885SU	Disposable Slightly Curved Metzenbaum Scissors Jaw Insert, 5 mm, 31 cm, 10 per box
1/1	P0886	Disposable Hook Scissors Jaw Insert, 5 mm, 31 cm, 10 per box
1/1	P0887	Disposable Mini-Metz Scissors, 5 mm, 31 cm, 10 per box
1/4	PO888 PO889	Disposable Metzenbaum Scissors Jaw Insert, 5 mm, 31 cm, 10 per box Disposable Metzenbaum Scissors Jaw Insert, 5 mm, 42 cm, 10 per box
10	PO958R	Quick-Snap™ Handle
	PM950R	Insulated Metal Handle

Monopolar Electrodes

Modular Monopolar Electrode System

QUALITY - German craftsmanship and design

DURABILITY - PEEK* insulation is resistant to cracks, peeling, shrinkage and malformations

MODULARITY - Exchange of tips and replacements as needed

ENHANCED SAFETY - Ceramic insulation to the working tip provides controlled energy delivery

SYSTEM COMPONENTS



Monopolar Handle - Complete instrument includes both Outer Tube and Inner Tube.

GK372R 5 mm diameter, 33 cm working length GK374R 5 mm diameter, 42 cm working length



Outer Tube

GK370P 5 mm, 33 cm **GK376P** 5 mm, 42 cm



Inner Tube and Handle

GK373R 5 mm, 33 cm **GK375R** 5 mm, 42 cm

MONOPOLAR EL	LECTRODE TIPS	S
--------------	---------------	---

-	GK383R Ceramic Electrode Tip, J-Hook
<u> </u>	GK384R Ceramic Electrode Tip, L-Hook
	GK385R Button Electrode Tip
	GK386R Spatula Electrode Tip

REUSABLE MONOPOLAR CABLES

	GK245	ERBE T Reusable Monopolar Cable, 4 mm diameter, 12 ft.
0	US354	VALLEYLAB, CONMED, BIRCHER, BOVIE, ERBE ICC Reusable Monopolar Cable, 4 mm diameter, 10 ft.

SPARE PARTS

PL170223 Contact Pin for HF Connection

GK299 Wire brush for cleaning the working tips

Suturing



INTRACORPOREAL		31 CM	37 CM	42 CM
1/1	Durogrip TC Needle Holder, straight, 5 mm	PL407R	PL404R	PL416R
1/1	Durogrip TC Needle Holder, curved left, 5 mm	PL408R	PL405R	PL417R
MIASS	Durogrip TC Needle Holder, curved right, 5 mm	PL413R	PL415R	
1/1	"The Viper" Durogrip TC Suture Assistant /Needle Holder, 2 x 3 teeth, 5 mm	PL409R		
1/1	Super Righting Needle Holder, 5 mm	PL414R		
1/1	Durogrip TC Needle Holder, Fenestrated, 5 mm	PL410R		
1/1	Durogrip TC Needle Holder, Crimped/Fenestrated, 5 mm	PL411R		
EXTRACORPOREAL				
1/1	CROCE/OLMI Dissecting and Suture Scissors with Knot Pusher, 5 mm	P0022R		
1/1	*Knot Pushing/Knot Tightening Forceps, 5 mm	PM654R PM973R PO958R		Jaw Insert Outer Tube Handle
	Knot Pusher, 5 mm, 33 cm	PL651R		
	Knot Pusher, 5 mm, 33 cm	PL432R		

^{*} Instrument sold by component

Bipolar Instruments



The enhanced safety of Ceramic Insulation in Reusable Bipolar Laparoscopic Modular Instruments

	COMPLETE INSTRUMENT	JAW INSERT	COMPONENTS INSULATED OUTER TUBES	STANDARD HANDLE
1/1	PM400R Bipolar Metzenbaum Scissors, 5 mm, 31 cm PM404R	PM430R	PM973R	PM450R
	Bipolar Metzenbaum Scissors, 5 mm, 42 cm	PM434R	PM981R	PM450R
1/1	PM407R Bipolar Mini Metzenbaum Scissors, 5 mm, 31 cm	PM437R	PM973R	PM450R
1/1	PM401R Bipolar Maryland Dissector, 5 mm, 31 cm	PM431R	PM973R	PM450R
all in	PM405R Bipolar Maryland Dissector, 5 mm, 42 cm	PM435R	PM981R	PM450R
1/1	PM408R Bipolar Maryland Dissector, fenestrated 5 mm, 31 cm PM415R	PM438R	PM973R	PM450R
	Bipolar Maryland Dissector, fenestrated 5 mm, 42 cm	PM445R	PM981R	PM450R
1/1	PM402R Bipolar Micro Forceps, 5 mm, 31 cm	PM432R	PM973R	PM450R
1/1	PM403R Bipolar Macro Forceps, 5 mm, 31 cm	PM433R	PM973R	PM450R
	PM406R Bipolar Macro Forceps, 5 mm, 42 cm	PM436R	PM981R	PM450R

Disposable Bipolar Instruments



		•
	ITEM NO.	DISPOSABLE BIPOLAR SCISSORS
	PL700	Metzenbaum Scissors, 5 mm, 31 cm, Six Single Sterile Scissors per Box
	PL701	Maryland Grasping Forceps, 5 mm, 31 cm, Six Single Sterile Graspers per Box
BIPOLAR CABLES & A	CCESSORIES	
/ ──	GN073	Bipolar Cable for Aesculap, Wolf, Martin, Berchtold, 12 ft.
	GK281	4 m length, suitable for units like: Aesculap, BOWA, EMC, Lamidey, Söring, Valleylab, 28.6 mm Pin Distance
	GK282	4 m length, suitable for units like: Codman CMC II, CMC III, 38 mm Pin Distance
	US349SF	Bipolar Cable for Valleylab, Conmed, Bovie, 12 ft. (Disposable, Single-use – 10 units/box)
GK216 Cable for Erbe		
GK217 Cable for Codman		

GK299 Wire Brush for cleaning jaws



	ITEM NO.	CLASSIC BIPOLAR TUBAL FORCEPS
-	MF970 MF980	32 cm, Syringe Handle Forceps 45 cm, Syringe Handle Forceps
	MF975 MF985	32 cm, Scissors Handle Forceps 45 cm, Scissors Handle Forceps
	MF971	32 cm, Jaw Insert Only
	MF981	45 cm, Jaw Insert Only
	MF966	Reusable Cable for Wolf Generator
	MF967	Reusable Cable for Valley Lab Generator

Specialty-Locking Graspers

Aesculap Slide-Lock Graspers

- Single- and double-action models
- Atraumatic and retraction jaw patterns
- Flush ports for optimal cleaning
- Excellent tissue grasping properties

ITEM NO.	DESCRIPTION	HANDLE	SHAFT
PL048R	Sovereign Atraumatic Grasper, single action, 5 mm, 42 cm	Slide Lock	Non-Rotatable
PLO49R	Sovereign Atraumatic Grasper, double action, 5 mm, 42 cm	Slide Lock	Rotatable
PL054R	Sovereign Atraumatic Grasper, single action, 5 mm, 36 cm	Slide Lock	Non-Rotatable
PL055R	Sovereign Atraumatic Grasper, double action, 5 mm, 36 cm	Slide Lock	Rotatable
PL056R	Sovereign Retraction Grasper, single action, 5 mm, 36 cm	Slide Lock	Non-Rotatable
PL057R	Sovereign Retraction Grasper, double action, 5 mm, 36 cm	Slide Lock	Rotatable
8360-00	Prestige Atraumatic Grasper, single action, 5 mm, 36 cm	Slide Lock	Non-Rotatable
8360-10	Prestige Atraumatic Grasper, double action, 5 mm, 36 cm	Slide Lock	Rotatable
8361-00	Prestige Retraction Grasper, single action, 5 mm, 36 cm	Slide Lock	Non-Rotatable
8361-10	Prestige Retraction Grasper, double action, 5 mm, 36 cm	Slide Lock	Rotatable



Sovereign® Graspers



Prestige[™] **Graspers**



Sovereign PL048R/54R Prestige 8360-00



Sovereign PL049R/55R Prestige 8360-10



Sovereign PL056R Prestige 8361-00



Sovereign PL057R Prestige 8361-10

Specialty Instruments Sets – Advanced Procedure

LAPAROSCOPIC VASCULAR CONTROL INSTRUMENTS



PL530R

Articulating Applicator for temporary endoscopic vessel clips 12.5 mm, 34 cm length


PL531R

Fixed Applicator for temporary endoscopic vessel clips, 12.5 mm, 34 cm length

PL664P

Rotation Knob, 5/package

Atraumatic Endo Vessel Clips (for temporary occlusion)

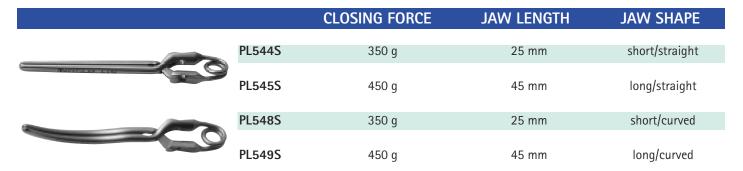

KIDNEY CLAMP

SIMON RENAL POLE CLAMP



Specialty Instruments Sets – Advanced Procedure

Optional Handles for Customizing Instruments – Interchangeable for 5 mm and 10 mm instruments









PO958R (Non-Ratchet)

PO959R (Ratchet)

PM950R (Non-Ratchet)

PM951R (Ratchet)







PM953R (In-Line)


PM954R (Ratchet)

PM957R (Ratchet) with lock

		COMPONENTS				
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES	
	PO014R Potts-DeMartell Scissors, 10 mm, 31 cm PM941R Protection Sleeve for safe	P0682R	PM992R	•	••	
	introduction through the 12.5 mm trocar					
1/2	PM171R* Satinsky Atraumatic Clamp, 38 mm jaw, 5 mm, 31 cm	PM711R	PM973R			
	PM177R* Satinsky Atraumatic Clamp, 51 mm jaw, 5 mm, 31 cm	PM616R	PM973R			
1/2	Order by component DeBakey Atraumatic Clamp, straight, 50 mm Jaw, 10 mm, 22 cm	PM712R	PM992R			
1/2	PM175R* DeBakey Atraumatic Clamp, curved, 5 mm, 31 cm	PM671R	PM973R		•	

^{*} For use with 12.5 mm flexible trocars

Optional Handles for Customizing Instruments – Interchangeable for 5 mm and 10 mm instruments

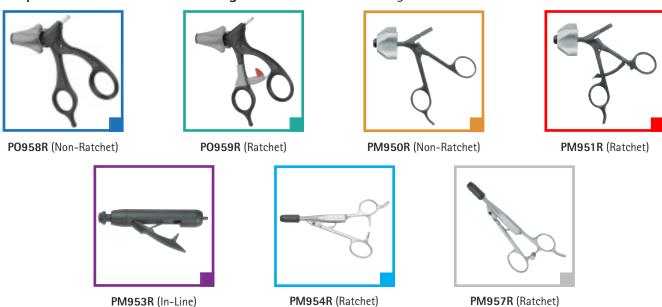

Instruments may be ordered COMPLETE or assembled from COMPONENTS ordered separately

with lock

		COMPONENTS				
	COMPLETE INSTRUMENT	JAW INSERT	INSULATED OUTER TUBES	STANDARD HANDLE	OPTIONAL HANDLES	
1/2	PM176R* DeBakey Atraumatic Clamp, S-shaped, 5 mm, 31 cm	PM689R	PM973R		•	
1/2	PM173R Crafoord Endo-vascular Clamp, short jaw, 5 mm, 31 cm PM178R Crafoord Endo-vascular Clamp, long jaw, 5 mm, 31 cm	PM713R PM719R	PM973R PM973R			
1/1	PO103R Overholt Dissecting Clamp, 90°, 10 mm, 31 cm	P0606R	PM975R			
1/1	PO218R Crile Dissecting Forceps, Curved 5 mm, 31 cm	P0771R	PM973R	•		
	PO213R Crile Dissecting Forceps, Curved 5 mm, 37 cm	P0766R	PM976R	•		

^{*} For use with 12.5 mm flexible trocars


Specialty Instruments Sets – Advanced Procedure


FLEXIBLE TROCARS



Complete trocar -100 mm length with valve and obturator


- Shaft sufficiently flexible to permit insertion of curved instruments and sufficiently firm to ensure access to the peritoneal orthoracic cavity
- For the use of curved instruments or special mouth shapes
- Trocar shaft in two different diameters and two different lengths

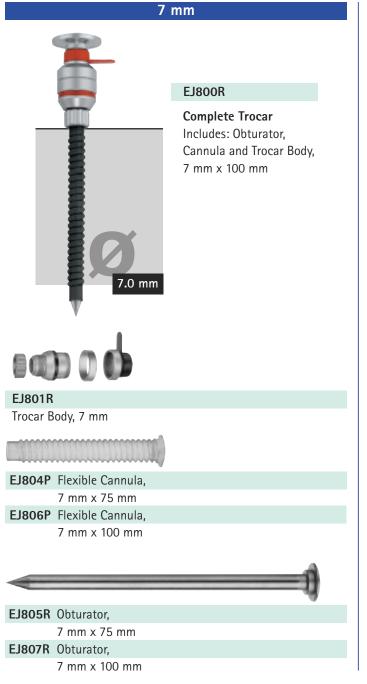

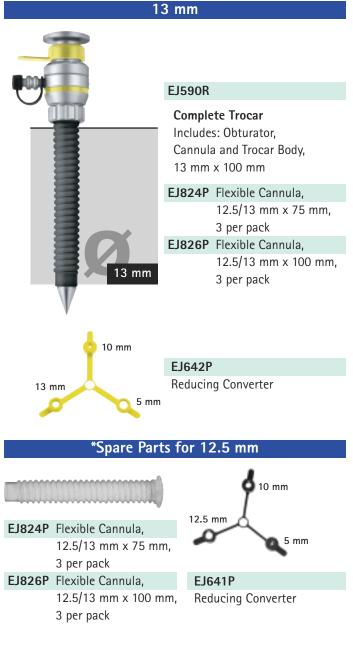



SLS Single-Fire Ligation System

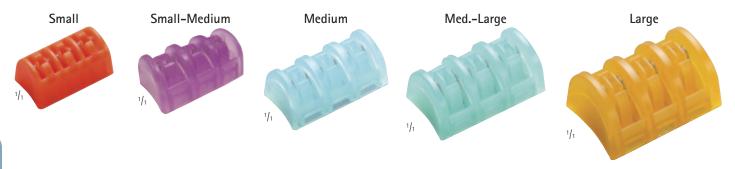

A Reliable and Cost-Efficient System that Delivers Superior Performance

- Clip cartridge is designed for easy loading and precise clip pickup
- Rotation knob allows for one finger, 360° rotation of applier shaft
- Modular clip applier can be disassembled for easy cleaning

ORDERING INFORMATION					
	Small	Small-Medium	Medium	Medium-Large	Large
Titanium Ligation Clips	J1180-1 Available in boxes of 30 cartridges of 6 clips = 180 clips	L5180-1 Available in boxes of 30 cartridges of 6 clips = 180 clips	B2180-1 Available in boxes of 30 cartridges of 6 clips = 180 clips	V3120-1 Available in boxes of 20 cartridges of 6 clips = 120 clips	O4120-1 Available in boxes of 20 cartridges of 6 clips = 120 clips
Laparoscopic Clip Appliers 13" (33 cm)	LCASLS1R Diameter 5.5 mm Length 30 cm	LCASLS2R Diameter 5.5 mm Length 30 cm	LCASLS3R Diameter 6.5 mm Length 30 cm LCASLS4R Diameter 10.0 mm Length 33 cm	LCASLS5R Diameter 10.0 mm Length 33 cm (shown left)	LCASLS6R Diameter 12.5 mm Length 33 cm

Challenger Ti-P Multi-Fire Ligation System

MULTI-FIRE

- Fast, precise and reliable clip application complimented by an ergonomic delivery
- Shafts available in 5 mm and 10 mm diameters
- Superior quality ligation clip available in S/M and M/L sizes
- Reusable clip applier, disposable clip cartridges
- Easy and simple handling

CHALLENGER TI-P PNEUMATIC CLIP APPLIERS



PL604R

- 5 mm diam., 31 cm length
 - 12 S/M clips per cartridge

PL606R

- 10 mm diam., 37 cm length
 8 M/L clips per cartridge
- 1/2

CLIP CARTRIDGES & CO2 CYLINDERS

S/M Clips

M/L Clips



- 12 sterile clip cartridges
 - 1 CO₂ cylinder per cartridge
 - For use with 5 mm clip applier

PL579T

- 12 sterile clip cartridges
 - 1 CO₂ cylinder per cartridge
 - For use with 10 mm clip applier





	Complete Shaft	PL522R		PL536R		PL538R	
		Shaft Length	Diameter	Shaft Length	Diameter	Shaft Length	Diameter
		310 mm	5 mm	370 mm	10 mm	260 mm	10 mm
	Handle	PL520R		PL520R		PL520R	
	Inner Tube	PL522801		PL536802		PL538802	
		Shaft Length		Shaft Length		Shaft Length	
		310 mm		370 mm		260 mm	
ţ	Push Rod	PL522800		PL536801		PL538801	
Parts		Shaft Length		Shaft Length		Shaft Length	
re		310 mm		370 mm		260 mm	
Spare	Locking Ring	PL520301		PL520301		PL520301	
•							
	CO ₂ Cylinder Only	PL575SU		PL575SU		PL575SU	
	Pack of 3						
	Single Sterile Packed						

Vascular Clip Applier/Remover



CLIP APPLIER/REMOVER



PL530R

Articulating Applicator for temporary endoscopic vessel clips, 12.5 mm, 34 cm length




PL531R

Fixed Applicator for temporary endoscopic vessel clips, 12.5 mm, 34 cm length

PL664P

Rotation Knob, 5/package

Atraumatic Endo Vessel Clips (for temporary occlusion)

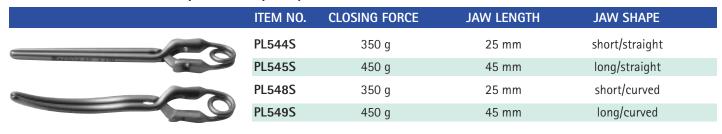

INTESTINAL BULLDOG CLIP



PL541S

Intestinal bulldog clip,

Jaw length: 70 mm, Closing force: 3.94 mm



HD Endoscopes



LAPAROSCOPES

ITEM NO.	DESCRIPTION	DIAMETER	VIEWING DIRECTION	LENGTH	Steam Autoclave	STERRAD®
PE030A	Needlescopic, Semi-rigid	2 mm	0°	30 cm	•	•
PE032A	Needlescopic, Semi-rigid	2.5 mm	0°	22 cm	•	•
PE485A	Full HD Endoscope	4 mm	0°	180 mm	•	•
PE505A	Full HD Endoscope	4 mm	30°	180 mm	•	•
PE522A	Full HD Endoscope	4 mm	0°	300 mm	•	•
PE508A	Full HD Endoscope	5 mm	30°	300 mm	•	•
PE590A	Full Screen Technology	5 mm	0°	30 cm	•	•
PE610A	Full Screen Technology	5 mm	30°	30 cm	•	•
PE889A	Full Screen Technology	10 mm	0°	33 cm	•	•
PE909A	Full Screen Technology	10 mm	30°	33 cm	•	•
PE969A	Full Screen Technology	10 mm	45°	33 cm	•	•
PE898A	Operative Laparoscope, with working channel for 5 mm instruments	10 mm	0°	27 cm	•	•

STORAGE SYSTEMS

JF431R	Single Scope Tray, 17", for 5 mm/10 mm diameter scopes
JF432R	Double Scope Tray, 17", for 5 mm/10 mm diameter scopes
JF433R	Double Scope Tray, 17", for 2 mm/3 mm diameter scopes
JF324R	Half-Size Storage Basket for 4 mm Perneczky angled scopes



Unitrac® Retraction & Holding System

The Unitrac is a pneumatically-powered retraction and scope holding system that allows you to experience a more efficient procedure and obstruction-free surgical site. Precise adjustments are possible with the push of a single button, thus freeing the surgical assistant to carry out higher priority tasks.

Features	Benefits
 Single-piece construction utilizes compressed air 	Easy setup and disassembly
source in the OR	
 Extensive range of retractor blade, hook, 	Applications in Minimally Invasive and Open
camera and endoscope holders	Surgical procedures
 Reposition with push of a single button 	Increases efficiency during procedures
 Maintains full holding force in event of OR 	Protects patient and OR personnel safety
compressed air failure	

UNITRAC ARM CONFIGURATIONS



OR Table Clamp





Compressed Air Hose



Y-connector for attaching 2 Unitrac arms to same compressed air source




Unitrac® Retraction & Holding System

The Unitrac is available in several configurations to ensure you have exactly what you need to achieve absolute rigid retraction and stable positioning. All configurations include the Unitrac pneumatic arm, a standard or flexible joint OR table clamp, and a compressed air hose or Y-connector for 2 Unitrac arms to be in use during the procedure.

Use the below ordering guide to select the Unitrac configuration that is right for you.

		ltem	No.	
	MB705R	MB706R	MB724R	MB725R
Unitrac Arm	•	•	•	•
OR Table Clamp	•	•		
OR Table Clamp w/ Flexible Joint			•	•
Compressed Air Hose	•		•	
Y-connector for attaching 2 Unitrac arms to same compressed air source		•		•

ENDOSCOPE HOLDERS		
	Item No.	Description
	RT046P	Adjustable holder for endoscopes with 3.0 - 7.5 mm diam., with PEEK insert RT055P
	RT055P	Replacement adjustable PEEK insert, 3.0 - 7.5 mm diam.
	RT081R	Holder for endoscopes with 10 mm diam., with PEEK insert RT086P
	RT086P	Replacement PEEK insert, 10 mm diam.

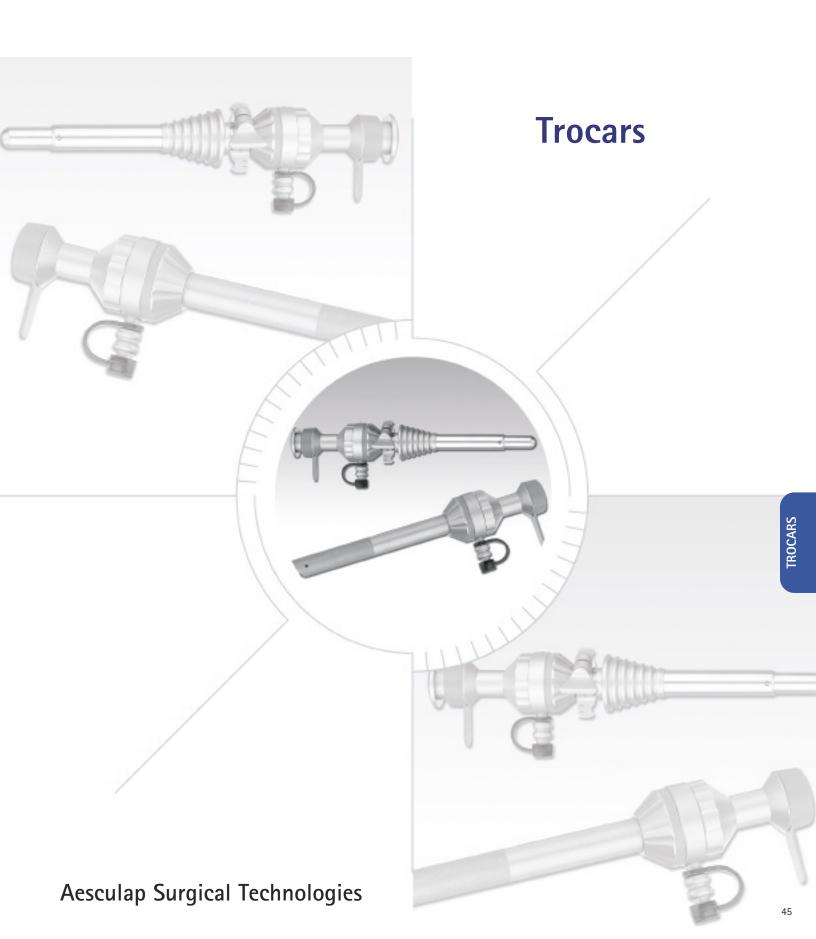

Reusable Trocar Systems


5.5 mm TROCARS





	WORKING LENGTH	COMPLETE TROCAR	CANNULA ONLY	PYRAMIDAL TIP OBTURATOR	CONICAL TIP OBTURATOR	BLUNT TIP OBTURATOR
400	60 mm		EJ702R	EJ705R	EJ706R	
	110 mm	EJ701R	EJ712R	EJ715R	EJ716R	EJ718R
Smooth Cannula	150 mm		EJ722R	EJ725R		EJ720R
	60 mm		EJ703R	EJ705R	EJ706R	
	110 mm		EJ713R	EJ715R	EJ716R	EJ718R
Threaded Cannula	150 mm					EJ720R
	60 mm		EJ704P	EJ705R	EJ706R	
	110 mm		EJ714P	EJ715R	EJ716R	EJ718R
Radiolucent Threaded Cannula	150 mm			EJ725R		EJ720R



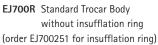




EJ710R Manual Opening Trocar Body

Sealing Caps for Reducing Sleeve 20/pkg


5.5 mm/10 mm	5.5 mm/12.5 mm	10 mm/12.5 mm
EJ601250	EJ603250	EJ605250

10 mm TROCARS

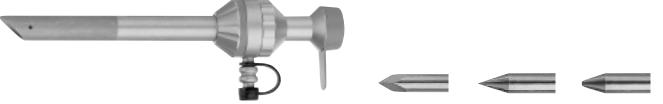

	WORKING LENGTH	COMPLETE TROCAR	CANNULA ONLY	PYRAMIDAL TIP OBTURATOR	CONICAL TIP OBTURATOR	BLUNT TIP OBTURATOR
	60 mm					
	110 mm	EJ751R	EJ752R	EJ755R	EJ756R	EJ558R
Smooth Cannula	150 mm		EJ762R	EJ765R	EJ766R	EJ568R
	60 mm					
	110 mm		EJ753R	EJ755R	EJ756R	EJ558R
Threaded Cannula	150 mm			EJ765R	EJ766R	EJ568R
	60 mm					
	110 mm		EJ754P	EJ755R	EJ756R	EJ558R
Radiolucent Threaded Cannula	150 mm			EJ765R	EJ766R	EJ568R



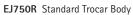



EJ760R Manual Opening Trocar Body



Reducing Sleeve

ITEM NO.	DIAMETER	LENGTH	GUIDE
F1601R	10 mm	180 mm	5.5 mm


Reducing Converter

ITEM NO.	TROCAR DIAMETER	CONVERTS TO:
EJ640P	10 mm & 12.5 mm	5 mm/10 mm

Sealing Caps for Reducing Sleeve 20/pkg

5.5 mm/10 mm	5.5 mm/12.5 mm	10 mm/12.5 mm
F1601250	F1603250	F1605250

Reusable Trocar Systems

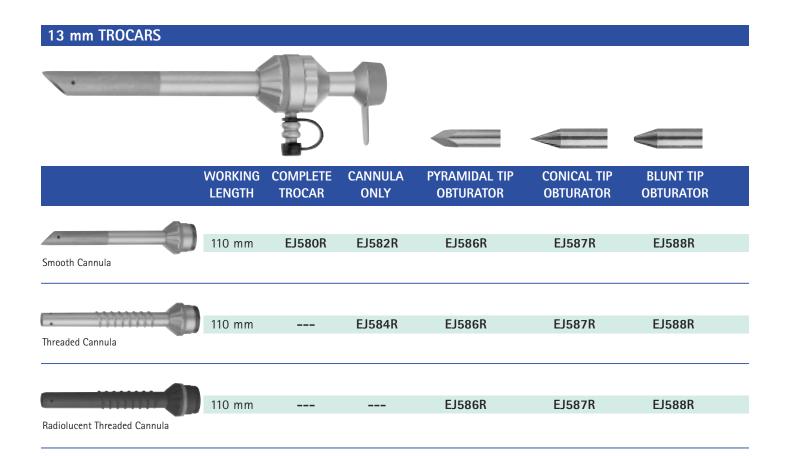



EJ581R Standard Trocar Body, without insufflation ring (order EJ750251 for insufflation ring)



Reducing Sleeve

ITEM NO.	DIAMETER	LENGTH	GUIDE
EJ602R	13 mm	180 mm	5.5 mm
EJ604R	13 mm	180 mm	10 mm



Reducing Converter

ITEM NO.	TROCAR DIAMETER	CONVERTS TO:
EJ642P	13 mm	5 mm, 10 mm, 13 mm

Sealing Caps for Reducing Sleeve 20/pkg

5.5 mm/10 mm	5.5 mm/12.5 mm	10 mm/12.5 mm	5.5 mm/13 mm	10 mm/13 mm
EJ601250	EJ603250	EJ605250	EJ646P	EJ645P

HASSON/MINI LAPAROTOMY STYLE TROCARS



COMPLETE TROCAR W/OBTURATOR

CANNULA ONLY

BLUNT TIP OBTURATOR



100 mm

EJ441R

EJ432R

EJ433R








EJ430R Stability Cone

TROCAR SPARE PARTS

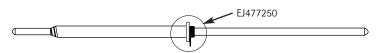

Dilation/Exchange Rod, 5 mm - 10 mm EJ477R EJ477250 Sealing Cap for EJ477R, 20/pkg.



Veress Needle

ITEM NO.	O.D. CANNULATED NEEDLE	O.D. CANNULA NEEDLE	LENGTH
PG003	1.6 mm	2.1 mm	100 mm
PG008	1.6 mm	2.1 mm	120 mm
PG011	1.6 mm	2.1 mm	150 mm

Reusable Trocar Systems

Trocar Spare Parts (continued)



EJ570P EJ571P EJ572P Silicone Flap Valve for 5.5 mm or 7 mm trocars, 20/pkg. Silicone Flap Valve for 10 mm or 12.5 mm trocars, 20/pkg.

Silicone Flap Valve for 13 mm trocars, 20/pkg.



EJ751251

Sealing Cap for Insufflation Connection for modular Trocars, 20/pkg.



EJ750251 EJ700251

Trocar Insufflation Ring with Luer Lock, 10 mm/12.5 mm, one piece Trocar Insufflation Ring with Luer Lock, 5 mm/7 mm, one piece



EJ750250 EJ700250 Trocar Ring without Luer Lock Connector, 10 mm/12.5 mm, one piece

Trocar Ring without Luer Lock Connector, 5 mm/7 mm, one piece



EJ702250

Plastic Ring for 5 mm trocars, 20/pkg. EJ752250 Plastic Ring for 10 mm trocars, 20/pkg.

EJ644P Plastic Ring for 13 mm trocars, 20/pkg.



SK271C

Metal Valve for Insufflation Connection, one piece



EJ446255 EJ417255 Sealing Cap for Modular Trocars, 5.5 mm, 20/pkg.


Sealing Cap for Modular Trocars, 7 mm, 20/pkg. EJ449255 Sealing Cap for Modular Trocars, 10 mm, 20/pkg.

EJ450255 Sealing Cap for Modular Trocars, 12.5 mm, 20/pkg. EJ643P Sealing Cap for Modular Trocars, 13 mm, 20/pkg.



MD156

Reducing Cap for Modular Trocars, 5/pkg.

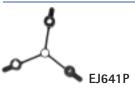

Reducing Converter; 5 mm, 10 mm, 12.5 mm


Flexible Trocars


FLEXIBLE TROCARS



Complete trocar -100 mm length with valve and obturator

- Shaft sufficiently flexible to permit insertion of curved instruments and sufficiently firm to ensure access to the peritoneal orthoracic cavity
- For the use of curved instruments or special mouth shapes
- Trocar shaft in two different diameters and two different lengths

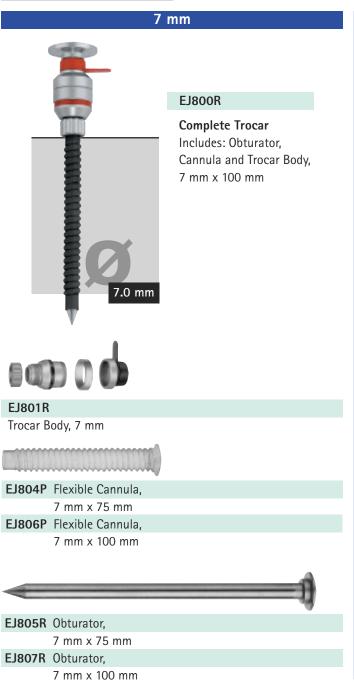

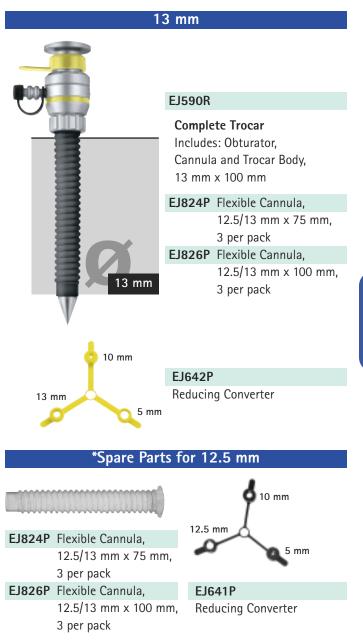




Light Sources

AXeL 180 XENON LIGHTSOURCE



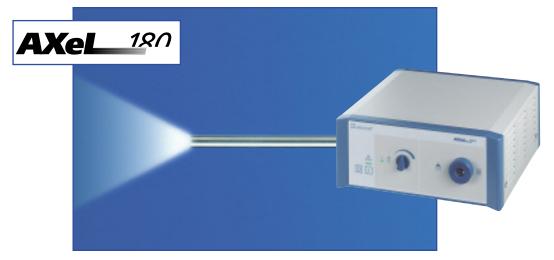

TECHNICAL DATA

LampXenonLamp Power180 WLamp Voltage15 V/12.0A

Light Intensity Continuously Adjustable

Color Temperature Approx. 6000 K

Dimensions (w x h x d) 12" x 5" x 12" (305 x 125 x 305 mm)

Weight 17 lbs., 17 oz. (6.5 kg)

ITEM NO. DESCRIPTION

OP930 AXeL Xenon Lightsource (180 W)

Consisting of:

OP931 Xenon Lamp, 15V, 180 W OP935 Standard Adapter, Aesculap

SYSTEM ACCESSORIES

OP936	Wolf LIGHT SOURCE Adapter
OP937	Olympus Lightsource Adapter
OP938	Circon ACMI Lightsource Adapter
OP913	Light Guide Cable, diameter: 5 mm, length: 250 cm (8.2 ft.) with additional bending protection
OP931	Xenon Spare Lamp

Light Cables

AESCULAP LIGHT CABLES and ADAPTERS



Light Cable Ordering Information

- Follow these steps to use the table below to identify the appropriate light cable:
 - Step 1 Select your light source from the horizontal column.
 - Step 2 Select your scope from the vertical column.
 - Step 3 The box where they meet represents the appropriate light cable.

(Example: Using an Aesculap light source and an Aesculap scope - Order MB111)



IF YOUR LIGHT SOURCE (LS) IS FROM:					
YOUR SCOPE (S) IS FROM	AESCULAP/ STORZ	OLYMPUS	CIRCON/ACMI	STRYKER	WOLF/DYONICS
AESCULAP/STORZ OLYMPUS	Light cable: MB111 • AESCULAP LS • AESCULAP S	Light cable: MB112 • OLYMPUS LS • AESCULAP S	Light cable: MB113 • CIRCON LS • AESCULAP S	Light cable: MB114 • STRYKER LS • AESCULAP S	Light cable: MB115 • WOLF LS • AESCULAP S
OLYMPUS (with condenser lens)	Light cable: MB121 • AESCULAP LS • OLYMPUS S	Light cable: MB122 • OLYMPUS LS • OLYMPUS S	Light cable: MB123 • CIRCON LS • OLYMPUS S	Light cable: MB124 • STRYKER LS • OLYMPUS S	Light cable: MB125 • WOLF LS • OLYMPUS S
CIRCON/ACMI/ STRYKER	Light cable: MB131 • AESCULAP LS • CIRCON S	Light cable: MB132 • OLYMPUS LS • CIRCON S	Light cable: MB133 • CIRCON LS • CIRCON S	Light cable: MB134 • STRYKER LS • CIRCON S	Light cable: MB135 • WOLF LS • CIRCON S
WOLF	Light cable: MB151 • AESCULAP LS • WOLF S	Light cable: MB152 • OLYMPUS LS • WOLF S	Light cable: MB153 • CIRCON LS • WOLF S	Light cable: MB154 • STRYKER LS • WOLF S	Light cable: MB155 • WOLF LS • WOLF S


• The universal light cable and all adapters can also be ordered separately

OP917

Light cable tester

MB100	Universal light cable with no end fittings, Light fiber bundle: 5 mm, Length: 10 ft.
MB001	Aesculap / Storz adapter for light cable (LS end)
MB002	Olympus adapter for light cable (LS end)
MB003	Circon / ACMI adapter for light cable (LS end)
MB004	Stryker adapter for light cable (LS end)
MB005	Wolf/Dyonics adapter for light cable (LS end)
MB010	Aesculap / Storz / Olympus adapter for light cable (S end)
MB030	Circon / ACMI / Stryker adapter for light cable (S end)
MB050	Wolf/Dyonics adapter for light cable (S end)



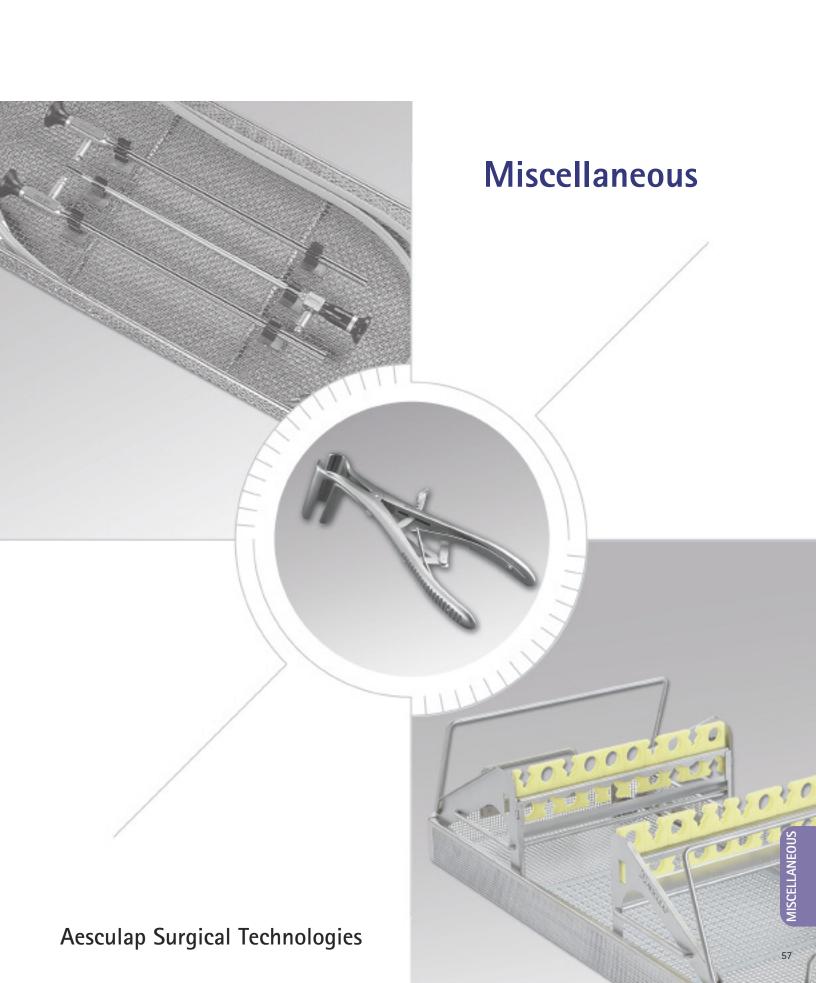

DISSECTING SPATULA

PL657R Dissecting Spatula, straight, 5 mm, 33 cm length

PL659R Dissecting Spatula, curved, 5 mm, 33 cm length (shown)

SCALPEL HANDLE for KNIFE BLADES Handle Features a retractable blade mechanism

BB060R Scalpel Handle, 5 mm, 33 cm length

KNIFE BLADES (order separately 10 per package)

.

BB363R

. .

BB364R

•

BB365R

BB367R

REUSABLE SUCTION IRRIGATION SYSTEMS



GF320R

Suction Tube

With Luer-Lock connection and finger valve. Can be used with Y-connector, SK330C, 5 mm diam., 33 cm length.

GF321R

Suction/Irrigation Tube

With Luer-Lock connection. Can be used with Y-connector, SK330C, 5 mm diam., 33 cm length.



SK330C

Y-Connector

For suction and irrigation in connection with SR550R, SR551R, SR552R, SR553R, GF320R and GF321R



Puncture Needle With Luer-Lock connection and finger valve, 5 mm diam., 33 cm length

SR551R

Injection Cannula, Sharp With Luer-Lock connection, 5 mm diam., 33 cm length

SR553R

Injection Cannula, Sharp With Luer-Lock connection, 5 mm diam., 45 cm length

LIGATURE INSTRUMENTS



PL650R

Loop ApplierOutside diam. – 5 mm
13 cm length
Inside diam. – 4.4 mm

PL650/250

Loop Applier Sealing Cap 20/pkg.



PL162R

Fascia Port-Closure Instrument 2.7 mm diam., 15.2 cm length, shaft only PL163R

ADDITIONAL ITEMS



PL430R

Myoma Fixation Screw 5 mm diam., 31 cm length



PL100R

Cholangiography Clamp with Catheter Guide

(Maximum catheter diameter-5 French), ratcheted, 5 mm diam., 33 cm length

PL401230

Sealing Cap for PL100R 20/pkg.



BV941R

Gall Bladder Extractor

Three-blade design, 55 mm blade length.

Use with 10 mm trocars.

Miscellaneous

ADDITIONAL ITEMS

BN320R

Exploratory Probe 5 mm diam., 40 cm length

BN321R

Exploratory Probe 5 mm diam., 50 cm length



PL406R

Anvil Grasper, 10 mm, 37 cm

LUBRICANTS



JG598

Aesculap Instrument Oil For box locks, ratchets, etc. Physiologically inert.



JG600

Aesculap Oil Spray
For maintenance of instruments prior to sterilization
Physiologically inert.

ENDOSCOPY TRAYS and BASKETS*



JF436R

Single Scope Tray, 10"
Two fit securely into
a Quarter Size container



JF432R or JF433R

Double Scope Tray, 17", with or without silicone supports Fits into a full size basket $17^{1}/2$ " long x $3^{3}/4$ " wide x $1^{3}/8$ " deep

ENDOSCOPY TRAYS and BASKETS*



MD322
Two Scope Basket
Includes: JF222R basket,
JG691 (4) instrument holding clamps

MD323


Three Scope Basket Includes: JF222R basket, JG691 (6) instrument holding clamps

MD324

Four Scope Basket Includes: JF222R basket, JG691 (8) instrument holding clamps



MD378
Endorack™ with Silicone
Cushioning Pad and Basket
Holds 14 endoscopic instruments

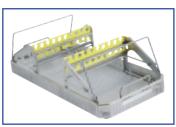

PL960R
Endorack, Full-Size
Storage System
Endorack with silicone mesh
mat and basket; holds 18
endoscopic instruments



JK444 - Solid Bottom
JN444 - Perforated Bottom
Full Size Container
Bottom with choice of lid:
 JK485 (red)
 JK486 (blue)
 JK487 (green)
 JK488 (gold)
 JK489 (silver)

^{*} For additional sterilization and storage products please contact Customer Service at 800-282-9000

IISCELLANEOUS

Miscellaneous

EXTRA LONG CONTAINER SYSTEM*




JN445 Perforated Bottom with Retention Plates



Lid with Retention Plates



MD379*
Endorack™ with Silicone
Cushioning Pad and Basket
Holds 14 endoscopic
instruments

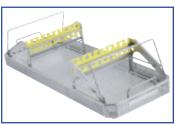

PL961R
Endorack, Extra-Long
Storage System
Endorack with silicone mesh
mat and basket; holds 18
endoscopic instruments

Product Index

Item No. Page No.	Item No. Page No.	Item No. Page No.
8360-0032	EJ587R48	EJ720R46
8360-1032	EJ588R48	EJ722R46
8361-0032	EJ590R36, 51	EJ725R46
8361-1032	EJ60125046, 47, 48	EJ75025050
B2180-138	EJ601R47	EJ75025148, 50
BB060R58	EJ602R48	EJ750R47, 49
BB363R58	EJ60325046, 47, 48	EJ75125150
BB364R58	EJ604R48	EJ751R47
BB365R58	EJ60525046, 47, 48	EJ75225050
BB367R58	EJ640P47	EJ752R47
BN320R60	EJ641P36, 50, 51	EJ753R47
BN321R60	EJ642P36, 48, 51	EJ754P47
BV941R59	EJ643P50	EJ755R47
EJ41725550	EJ644P50	EJ756R47
EJ430R49	EJ645P48	EJ760R47
EJ432R49	EJ646P48	EJ762R47
EJ433R49	EJ70025050	EJ765R47
EJ441R49	EJ70025146, 50	EJ766R47
EJ44625550	EJ700R46	EJ800R36, 51
EJ44925550	EJ701R46	EJ801R36, 51
EJ45025550	EJ70225050	EJ804P36, 51
EJ47725049	EJ702R46	EJ805R36, 51
EJ477R49	EJ703R46	EJ806P36, 51
EJ558R47	EJ704P46	EJ807R36, 51
EJ568R47	EJ705R46	EJ824P36, 51
EJ570P50	EJ706R46	EJ826P36, 51
EJ571P50	EJ710R46	GF320R58
EJ572P50	EJ712R46	GF321R58
EJ580R48	EJ713R46	GK21631
EJ581R48	EJ714P46	GK21731
EJ582R48	EJ715R46	GK24526, 28
EJ584R48	EJ716R46	GK24626
EJ586R48	EJ718R46	GK28131

Item No. Page No.	Item No. Page No.	Item No. Page No.
GK28231	L5180-138	MB15355
GK29928, 31	LCASLS1R38	MB15455
GK370P28	LCASLS2R38	MB15555
GK372R28	LCASLS3R38	MB705R44
GK373R28	LCASLS4R38	MB706R44
GK374R28	LCASLS5R38	MB724R44
GK375R28	LCASLS6R38	MB725R44
GK376P28	MB00155	MD15650
GK383R28	MB00255	MD32261
GK384R28	MB00355	MD32361
GK385R28	MB00455	MD32461
GK386R28	MB00555	MD37861
GK469R26	MB01055	MD37962
GN07331	MB03055	MF96631
J1180-138	MB05055	MF96731
JF222R61	MB10055	MF97031
JF324R42	MB11155	MF97131
JF431R42	MB11255	MF97531
JF432R42, 60	MB11355	MF98031
JF433R42, 60	MB11455	MF98131
JF436R60	MB11555	MF98531
JG59860	MB12155	04120-138
JG60060	MB12255	OP91354
JG69161	MB12355	OP91755
JK44461	MB12455	0P93054
JK48561	MB12555	OP93154
JK48661	MB13155	OP93554
JK48761	MB13255	OP93654
JK48861	MB13355	OP93755
JK48961	MB13455	OP93854
JK49062	MB13555	PE030A42
JN44461	MB15155	PE032A42
JN44562	MB15255	PE485A42

Product Index

PE505A 42 PL413R 29 PL657R 58 PE508A 42 PL414R 29 PL659R 58 PE52A 42 PL416R 29 PL66AP 33 40 PE590A 42 PL417R 29 PL700 31 7 PE68B/219 42,55 PL430R 59 PL960R 61 PE68B/221 42,55 PL432R 29 PL961R 62 PE889A 42 PL52301 39 PM015R 12 PE898A 42 PL522800 39 PM16R 15 PE969A 42 PL522800 39 PM178R 35 PE969A 42 PL52280 39 PM178R 35 P6003 49 PL52280 39 PM178R 35 P6004 42 PL530R 33,40 PM178R 35 P6008 49 PL530R 33,40 PM176R 35 PL	Item No. Page No.	Item No. Page No.	Item No. Page No.
PE522A	PE505A42	PL413R29	PL657R58
PE590A	PE508A42	PL414R29	PL659R58
PE610A 42 PL417R 29 PL701 31 PE688/219 42,55 PL430R 59 PL960R 61 PE688/221 42,55 PL432R 29 PL961R 62 PE889A 42 PL520301 39 PM015R 12 PE898A 42 PL52080 39 PM16SR 13 PE909A 42 PL522800 39 PM16SR 15 PE969A 42 PL522800 39 PM173R 34 PG003 49 PL522R 39 PM173R 35 PG008 49 PL530R 33, 40 PM175R 34 PG011 49 PL531R 33, 40 PM176R 35 PL048R 32 PL536801 39 PM177R 34 PL049R 32 PL536801 39 PM407R 35 PL051R 32 PL536801 39 PM407R 30 PL057R 32	PE522A42	PL415R29	PL664P33, 40
PE688/21942, 55 PL430R59 PL960R61 PE688/22142, 55 PL432R29 PL961R62 PE889A42 PL52030139 PM015R12 PE898A42 PL520R39 PM16R13 PE99A42 PL52280039 PM16SR15 PE969A42 PL52280139 PM171R34 PG00349 PL532R39 PM173R35 PG00849 PL53R33,40 PM175R34 PG01149 PL531R33,40 PM176R35 PL048R32 PL53680139 PM177R34 PL049R32 PL53680139 PM400R35 PL055R32 PL53680139 PM400R30 PL056R32 PL53880139 PM401R30 PL057R32 PL53880239 PM402R30 PL057R32 PL548S33,40 PM404R30 PL100R59 PL541S40 PM404R30 PL102232 PL548S33,40 PM405R30 PL1702332 PL548S33,40 PM407R30 PL17025326 PL548S33,40 PM406R30 PL404R29	PE590A42	PL416R29	PL70031
PE688/221 42,55 PL432R 29 PL961R 62 PE889A 42 PL520301 39 PM015R 12 PE898A 42 PL52080 39 PM016R 13 PE909A 42 PL522801 39 PM178R 34 PG003 49 PL522R 39 PM173R 35 PG008 49 PL530R 33, 40 PM175R 34 PG011 49 PL531R 33, 40 PM176R 35 PL048R 32 PL536801 39 PM177R 34 PL049R 32 PL536802 39 PM178R 35 PL049R 32 PL536803 39 PM400R 30 PL05R 32 PL538803 39 PM401R 30 PL05FR 32 PL538803 39 PM401R 30 PL05FR 32 PL538803 39 PM402R 30 PL162F 32	PE610A42	PL417R29	PL70131
PE889A	PE688/21942, 55	PL430R59	PL960R61
PE898A	PE688/22142, 55	PL432R29	PL961R62
PE909A	PE889A42	PL52030139	PM015R12
PE969A	PE898A42	PL520R39	PM016R13
P6003	PE909A42	PL52280039	PM165R15
PG008 49 PL530R 33, 40 PM175R 34 PG011	PE969A42	PL52280139	PM171R34
PG011 49 PL531R	PG00349	PL522R39	PM173R35
PL048R 32 PL536801 39 PM177R 34 PL049R 32 PL536802 39 PM178R 35 PL054R 32 PL536R 39 PM400R 30 PL055R 32 PL538801 39 PM401R 30 PL057R 32 PL538R 39 PM402R 30 PL100R 59 PL541S 40 PM404R 30 PL162R 59 PL544S 33, 40 PM405R 30 PL170223 28 PL545S 33, 40 PM406R 30 PL170251 26 PL548S 33, 40 PM407R 30 PL401230 59 PL549S 33, 40 PM408R 30 PL404R 29 PL575SU 39 PM408R 30 PL405R 29 PL575SU 39 PM430R 30 PL406R 60 PL579T 39 PM431R 30 PL408R 29 PL606R 39 PM433R 30 PL409R 29 <td< td=""><td>PG00849</td><td>PL530R33, 40</td><td>PM175R34</td></td<>	PG00849	PL530R33, 40	PM175R34
PL049R 32 PL536802 39 PM400R 30 PL054R 32 PL538801 39 PM401R 30 PL056R 32 PL538802 39 PM402R 30 PL057R 32 PL538R 39 PM402R 30 PL106R 59 PL541S 40 PM404R 30 PL162R 59 PL544S 33, 40 PM405R 30 PL170223 28 PL545S 33, 40 PM406R 30 PL170251 26 PL548S 33, 40 PM407R 30 PL401230 59 PL549S 33, 40 PM408R 30 PL404R 29 PL575SU 39 PM408R 30 PL405R 29 PL575SU 39 PM431R 30 PL406R 60 PL579T 39 PM431R 30 PL407R 29 PL606R 39 PM432R 30 PL409R 29 PL606R	PG01149	PL531R33, 40	PM176R35
PL054R 32 PL536R 39 PM400R 30 PL055R 32 PL538801 39 PM401R 30 PL056R 32 PL538802 39 PM402R 30 PL057R 32 PL538R 39 PM403R 30 PL100R 59 PL541S 40 PM404R 30 PL162R 59 PL544S 33, 40 PM405R 30 PL17023 28 PL548S 33, 40 PM406R 30 PL170251 26 PL548S 33, 40 PM407R 30 PL401230 59 PL549S 33, 40 PM408R 30 PL404R 29 PL574T 39 PM408R 30 PL405R 29 PL575SU 39 PM430R 30 PL406R 60 PL579T 39 PM431R 30 PL407R 29 PL604R 39 PM432R 30 PL409R 29 PL606R 39 PM434201 26 PL410R 29 P	PL048R32	PL53680139	PM177R34
PL055R	PL049R32	PL53680239	PM178R35
PL056R	PL054R32	PL536R39	PM400R30
PL057R 32 PL538R 39 PM403R 30 PL100R 59 PL541S 40 PM404R 30 PL162R 59 PL544S 33, 40 PM405R	PL055R32	PL53880139	PM401R30
PL100R 59 PL541S 40 PM404R 30 PL162R 59 PL544S 33, 40 PM405R 30 PL170223 28 PL545S 33, 40 PM406R 30 PL170251 26 PL548S 33, 40 PM407R 30 PL401230 59 PL549S 33, 40 PM408R 30 PL404R 29 PL574T 39 PM415R 30 PL405R 29 PL575SU 39 PM430R 30 PL407R 29 PL604R 39 PM431R 30 PL408R 29 PL606R 39 PM433R 30 PL409R 29 PL650/250 59 PM434201 26 PL410R 29 PL650R 59 PM434R 30	PL056R32	PL53880239	PM402R30
PL162R59 PL544S33, 40 PM405R30 PL17022328 PL545S33, 40 PM406R30 PL17025126 PL548S33, 40 PM407R30 PL40123059 PL549S33, 40 PM408R30 PL404R29 PL574T39 PM415R30 PL405R29 PL575SU39 PM430R30 PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R	PL057R32	PL538R39	PM403R30
PL17022328 PL545S33, 40 PM406R30 PL17025126 PL548S33, 40 PM407R30 PL40123059 PL549S33, 40 PM408R30 PL404R29 PL574T39 PM415R30 PL405R29 PL575SU39 PM430R30 PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL100R59	PL541S40	PM404R30
PL17025126 PL548S33, 40 PM407R30 PL40123059 PL549S33, 40 PM408R30 PL404R29 PL574T39 PM415R30 PL405R29 PL575SU39 PM430R30 PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL162R59	PL544S33, 40	PM405R30
PL40123059 PL549S33, 40 PM408R30 PL404R29 PL574T39 PM415R30 PL405R29 PL575SU39 PM430R30 PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL17022328	PL545S33, 40	PM406R30
PL404R	PL17025126	PL548S33, 40	PM407R30
PL405R29 PL575SU39 PM430R30 PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL40123059	PL549S33, 40	PM408R30
PL406R60 PL579T39 PM431R30 PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL404R29	PL574T39	PM415R30
PL407R29 PL604R39 PM432R30 PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL405R29	PL575SU39	PM430R30
PL408R29 PL606R39 PM433R30 PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL406R60	PL579T39	PM431R30
PL409R29 PL650/25059 PM43420126 PL410R29 PL650R59 PM434R30	PL407R29	PL604R39	PM432R30
PL410R29 PL650R59 PM434R30	PL408R29	PL606R39	PM433R30
	PL409R29	PL650/25059	PM43420126
PL411R29 PL651R29 PM435R30	PL410R29	PL650R59	PM434R30
	PL411R29	PL651R29	PM435R30

Item No. Page No.	Item No. Page No.	Item No. Page No.
PM436R30	PM953R12, 13, 14, 15, 16, 17, 18,	PO103R16, 35
PM437R30	19, 20, 21, 22, 23, 34, 35	PO104R19
PM438R30	PM954R34, 35	PO105R20
PM445R30	PM957R34, 35	PO107R17
PM450R30	PM973R12, 13, 14, 15, 16, 17, 18,	PO109R17
PM60020126	19, 20, 21, 22, 23, 26, 29,	PO111R19
PM604R22	30, 34, 35	PO112R17
PM609R22	PM974R18, 19, 21, 22, 23, 26	PO113R18
PM612R23	PM975R16, 20, 26, 35	PO114R19
PM614R22	PM976R15, 17, 26, 35	PO117R19
PM616R34	PM977R18, 19, 23, 26	PO118R18
PM61820126	PM978R26	PO128R19
PM621R22	PM981R12, 13, 14, 15, 17, 18, 19,	P0129R14
PM654R29	20, 21, 26, 30	P0130R20
PM657R23	PM982R19, 26	P0132R20
PM671R34	PM992R34	P0134R20
PM679R23	PM995R26	P0137R18
PM688R23	P0001R12	PO141R16
PM689R35	P0002R13	P0147R16
PM697R12	P0003R12	P0149R15
PM698R13	P0004R12	P0151R16
PM699R15	P0006R13	P0187R16
PM711R34	P0011R13	P0190R17
PM712R34	P0014R34	PO191R18
PM713R35	P0022R14, 29	PO196R19
PM719R35	P0025R12	P0199R14
PM941R34	P0027R12	P0200R14
PM95025026	P0028R13	P0201R20
PM950R12, 13, 14, 15, 16, 17, 18,	P0090R13	P0202R15
19, 20, 21, 22, 23, 24, 27,	P0091R13	P0203R17
34, 35	P0100R15	P0204R17
PM951R12, 13, 14, 15, 16, 17, 18,	PO101R14	P0205R19
19, 20, 21, 22, 23, 24, 34, 35	PO102R15	PO206R15

Product Index

P0207R14 P0602R12 P0682R34 P0208R15 P0603R12 P0705R21 P0209R18 P0605R15 P0727R14 P0210R14 P0606R16, 35 P0730R16 P0211R15 P0607R14 P0735R17 P0212R	
P0209R18 P0605R15 P0727R14 P0210R14 P0606R16, 35 P0730R16 P0211R15 P0607R14 P0735R17 P0212R18 P0608R15 P0736R18 P0213R15, 35 P0610R19 P0741R19 P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0210R14 P0606R16, 35 P0730R16 P0211R15 P0607R14 P0735R17 P0212R18 P0608R15 P0736R18 P0213R15, 35 P0610R19 P0741R19 P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0211R15 P0607R14 P0735R17 P0212R18 P0608R15 P0736R18 P0213R15, 35 P0610R19 P0741R19 P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0212R18 P0608R15 P0736R18 P0213R15, 35 P0610R19 P0741R19 P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0213R15, 35 P0610R19 P0741R19 P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0214R17 P0611R20 P0744R12 P0215R18 P0613R17 P0745R12	
P0215R18 P0745R12	
P0216R17 P0615R17 P0746R12	
P0217R17 P0617R19 P0747R14	
P0218R15, 35 P0618R17 P0748R14	
P0221R20 P0749R20	
P0222R20 P0620R19 P0750R21	
P0223R20 P0751R21	
P0224R21 P0624R18 P0752R15	
P0225R21 P0753R17	
P0228R33 P0626R21 P0754R13	
P0229R19 P0755R17	
P0230R20 P0628R14 P0756R21	
P0231R16 P0757R19	
P0500R21 P0646R20 P0758R15	
P0501R21 P0647R19 P0759R21	
P0502R21 P0648R20 P0760R14	
P0505R20 P0761R15	
P0509R21 P0762R18	
P0510R21 P0763R14	
P0511R21 P0657R18 P0764R15	
P0512R16 P0765R18	
P0550R14 P0667R16 P0766R15, 35	
P0559R14 P0767R17	
P0600R12 P0678R16 P0768R18	
P0601R13 P0769R17	

Item No.	Page No.	Item No.	Page No.
P0770R	17	V3120-1	38
P0771R	15, 35		
P0774R	20		
P0775R	20		
P0776R	20		
P0777R	21		
P0778R	21		
P0782R	20		
P0783R	19		
P0784R	16		
P0790R	13		
P0791R	13		
P0885SU	27		
P0886	27		
P0887	27		
P0888	27		
P0889	27		
P0950400	26		
P0958R	12, 13, 14, 15, 16, 17, 18,		
	19, 20, 21, 22, 23, 25, 27,		
	29, 34, 35		
P0959R	12, 13, 14, 15, 16, 17, 18,		
	19, 20, 21, 22, 23, 25, 34, 35		
RT046P	44		
RT055P	44		
SK271C	50		
SK330C	58		
SR550R	58		
SR551R	58		
SR552R	58		
SR553R	58		
US349SP	31		
US354	26, 28		

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogs and price lists and to take legal actions. ©2010 AESCULAP. ALL RIGHTS RESERVED. PRINTED IN THE USA. Aesculap is an equal opportunity employer	