

Top



UNIVERSAL 320/320 R
Benchtop Centrifuges
classic/cooled

There are good reasons why we have called these centrifuges UNIVERSAL. Their excellent performance and comprehensive range of accessories enable them to carry out virtually any centrifuging task. In the low speed range, the UNIVERSAL 320 and 320 R process standard tubes up to a volume of 100 ml, microtitre plates, conical tubes with screw cap, blood collection tubes and cyto accessories. Microlitre tubes can be accelerated to an RCF of 21,382. Centrifuge tubes taking higher volumes (up to 85 ml) also profit from the extraordinary motor power of the UNIVERSAL 320 / 320 R. At an RCF of 9,509, bacteria and yeast sediment quickly and easily.



UNIVERSAL 320 classic
Cat. No. 1401

PERFORMANCE

- High RCF
 - 85 ml tubes up to 9,509
 - 15 ml tubes up to 16,582
 - Microlitre tubes up to 21,382
- Extremely short run-up and run-down times

DESIGN

- Metal housing
- Metal lid
- Stainless steel centrifuging chamber
- Viewing port in the lid

USER-FRIENDLINESS

- Powered lid-locking
- Quick-entry foil keypad
- Easily exchangeable rotors
- 10 programmable memories

REFRIGERATION (UNIVERSAL 320 R)

- Infinitely variable setting from – 20 °C to + 40 °C
- Fast cool function

SAFETY

- Lid locking and holding during rotor run
- Emergency lid-lock release
- Motor-overheating protection
- Chamber-overheating protection (with the UNIVERSAL 320 R)
- Imbalance switch-off
- Automatic rotor recognition
- Lid dropping protection

OPTIONS

As options, the UNIVERSAL 320 R is also available

- with connection to an external refrigerated/heating circulator (provided by the customer) – **Cat. No. 1406-50** (without internal cooling unit). Depending on circulator performance, the centrifuge can heat or cool within a range of 0 °C to 90 °C,
- with connection for flushing with nitrogen – **Cat. No. 1406-20** (without monitoring and controlling system).

The display panel and all controls on the UNIVERSAL 320 and 320 R centrifuges have been designed ergonomically. Information on the display is clearly positioned and perfectly legible. During centrifugation, the actual parameter values are indicated. Nine acceleration and braking ramps, or unbraked run-down, can be pre-selected. Ramp 9 always corresponds to the shortest possible run-up or run-down time. The stored parameter combination remains in the memory even after the centrifuge has been switched off. By entering the rotor radius r/mm, the centrifuge calculates and indicates the rotational speed RPM or the relative centrifugal force RCF. The twist knob allows quick and easy entry or changing of the parameters.

DIGITAL DISPLAY AND CONTROL PANEL

KEYPAD

	Precooling (UNIVERSAL 320 R)
RCF	Change-over from RPM to RCF indication
SELECT	Selection of parameters
START IMPULS	Starts centrifugation / Short cycle mode
STOP OPEN	Stops centrifugation manually / Opens the lid

ENTRY OF THE PARAMETERS

P	Entry of the program number, 10 programmable memories
T/°C	Entry of the temperature in increments of 1 within a range of - 20 °C to + 40 °C (UNIVERSAL 320 R)
RCF	Entry in increments of 10
RPM	Entry in increments of 10
r/mm	Entry of the radius in mm
t/min	Entry of the centrifuging time in minutes and seconds (max. 99 min : 59 sec)
	Entry of the acceleration ramp 1-9 Entry of the braking ramp 1-9



UNIVERSAL 320 R control panel
(N Plus control panel)



UNIVERSAL 320 R cooled
Cat. No. 1406

ROTORS AND ACCESSORIES

UNIVERSAL 320 / 320 R centrifuges accommodate rotors and accessories for virtually every centrifuging task required in industrial, research and clinical laboratories.

Quick overview of rotors and suitable accessories. The figures highlighted under the arrows refer to the pages where technical and ordering details are listed.

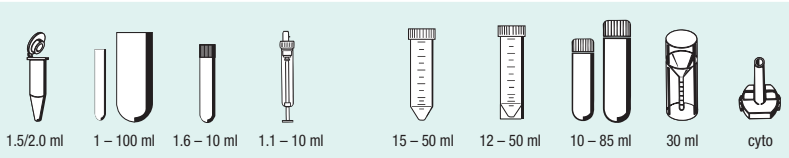
Tubes / plates	Max. capacity
standard tubes	32 x 15 ml or 4 x 100 ml
hematocrit capillaries	24
conical tubes with screw cap	32 x 15 ml or 8 x 50 ml
microtitre plates	10
filter plates	2
cyto suspensions	8
blood collection tubes	48
round-bottom tubes with screw cap	6 x 85 ml
microlitre tubes	30 x 1.5 / 2.0 ml

Rotors	Max. RCF
swing-out rotor	4,193
swing-out rotor for conical tubes with screw cap	3,857
swing-out rotor for microtitre and culture plates	2,218
swing-out rotor for 8 cyto suspensions	2,218
angle rotor for round-bottom tubes with screw cap (6 x 85 ml)	9,509
angle rotors (12-, 24- or 30-place) for microlitre tubes (0.2 – 2.0 ml)	21,382



Swing-out rotor 4-place
n = 5,000 min⁻¹
max. RCF 4,193
Cat. No. 1494

some accessories with biocontainment⁶⁾



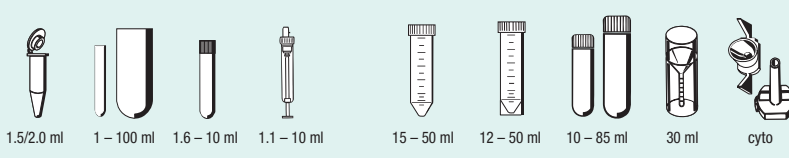
1.5/2.0 ml 1 – 100 ml 1.6 – 10 ml 1.1 – 10 ml 15 – 50 ml 12 – 50 ml 10 – 85 ml 30 ml cyto

→ Page 6-8



Swing-out rotor 4-place
n = 4,000 min⁻¹
max. RCF 2,665
Cat. No. 1624

some accessories with biocontainment⁶⁾



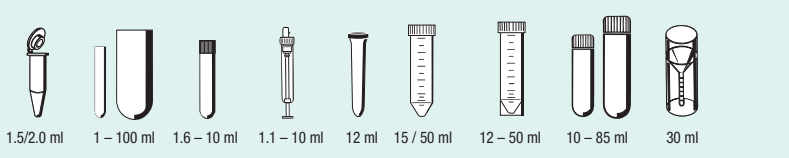
1.5/2.0 ml 1 – 100 ml 1.6 – 10 ml 1.1 – 10 ml 15 – 50 ml 12 – 50 ml 10 – 85 ml 30 ml cyto

→ Page 9-10



Swing-out rotor 4-place
n = 4,500 min⁻¹
max. RCF 3,328
Cat. No. 1324

some accessories with biocontainment⁶⁾

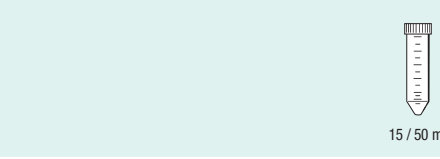


1.5/2.0 ml 1 – 100 ml 1.6 – 10 ml 1.1 – 10 ml 12 ml 15 / 50 ml 12 – 50 ml 10 – 85 ml 30 ml

→ Page 11-13



Swing-out rotor 6-place
n = 4,000 min⁻¹
max. RCF 2,701
Cat. No. 1619

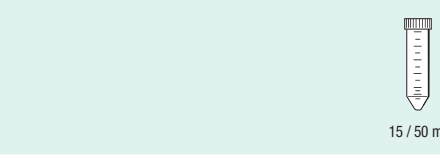


15 / 50 ml

→ Page 13

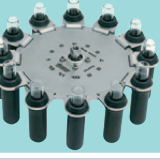


Swing-out rotor 8-place
n = 5,000 min⁻¹
max. RCF 3,857
Cat. No. 1617

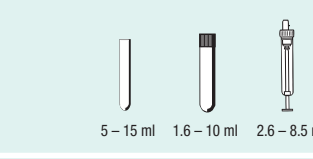


15 / 50 ml

→ Page 13



Swing-out rotor 12-place
n = 5,000 min⁻¹
max. RCF 4,193
Cat. No. 1628

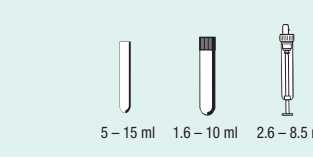


5 – 15 ml 1.6 – 10 ml 2.6 – 8.5 ml

→ Page 13



Swing-out rotor 8-place
n = 4,000 min⁻¹
max. RCF 2,415
Cat. No. 1611

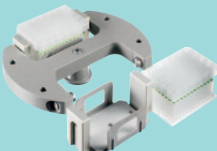
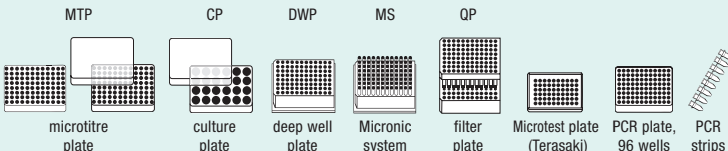
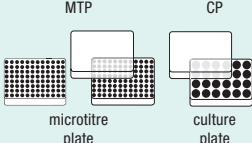
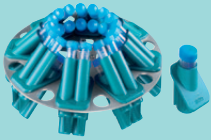
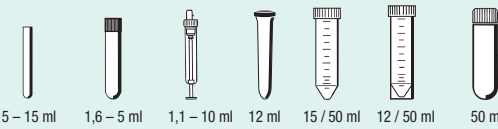
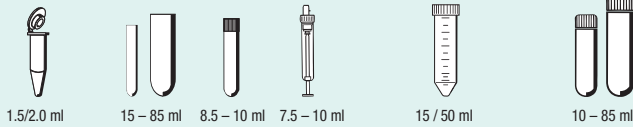
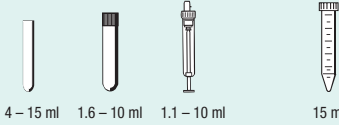
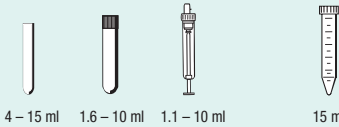
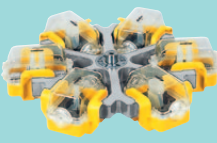


5 – 15 ml 1.6 – 10 ml 2.6 – 8.5 ml

→ Page 14

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

ROTORS AND ACCESSORIES

 Swing-out rotor 2-place ↗90° n = 4,000 min⁻¹ max. RCF 2,218 Cat. No. 1460	 MTP CP DWP MS QP microtitre plate culture plate deep well plate Micronic system filter plate Microtest plate (Terasaki) PCR plate, 96 wells PCR strips → Page 14
 Swing-out rotor 4-place ↗90° n = 4,000 min⁻¹ max. RCF 1,896 Cat. No. 1645	 MTP CP microtitre plate culture plate → Page 15
 Angle rotor 12-place ↗45° n = 15,000 min⁻¹ max. RCF 16,602 Cat. No. 1612	 0.2 – 2.0 ml → Page 15
auto-clavable optional with biocontainment⁶⁾  Angle rotor 24-place ↗40° n = 15,000 min⁻¹ max. RCF 21,382 Cat. No. 1420-B	 0.2 – 2.0 ml → Page 15
auto-clavable with bio-containment⁶⁾  Angle rotor 30-place ↗45° n = 14,000 min⁻¹ RCF 21,255 Cat. No. 1689-A	 0.2 – 2.0 ml → Page 16
 Angle rotor 8-place ↗45° n = 4,500 min⁻¹ max. RCF 3,283 Cat. No. 1418	 5 – 15 ml 1,6 – 5 ml 1,1 – 10 ml 12 ml 15 / 50 ml 12 / 50 ml 50 ml → Page 16
 Angle rotor 6-place ↗35° n = 9,000 min⁻¹ max. RCF 9,509 Cat. No. 1620A	 1.5/2.0 ml 15 – 85 ml 8.5 – 10 ml 7.5 – 10 ml 15 / 50 ml 10 – 85 ml → Page 17
 Angle rotor 12-place ↗35° n = 6,000 min⁻¹ max. RCF 4,146 Cat. No. 1613	 4 – 15 ml 1.6 – 10 ml 1.1 – 10 ml 15 ml → Page 17
 Angle rotor 12-place ↗35° n = 12,000 min⁻¹ max. RCF 16,582 Cat. No. 1615	 4 – 15 ml 1.6 – 10 ml 1.1 – 10 ml 15 ml → Page 18
 Hematocrit rotor 24-place n = 15,000 min⁻¹ max. RCF 21,382 Cat. No. 1650	 capillary → Page 18
 Cyto rotor 6-place n = 4,000 min⁻¹ max. RCF 2,039 Cat. No. 1626	Detailed information on our comprehensive range of cyto accessories is provided in a separate brochure. We will be pleased to send one to you on request.  cyto → Page 19
 Cyto rotor 8-place n = 4,000 min⁻¹ max. RCF 2,218 Cat. No. 1648	 cyto → Page 19

Additional tubes e.g. stool tubes or spin columns on request.

ROTORS AND ACCESSORIES

Swing-out rotor, 4-place



illustrated with carriers 1425

Cat. No. (without carriers) **1494**

capacity in ml	1.5	2.0	1	3	5	6	7	9	15	25	50	100	suspension for cyto accessories
Ø x L in mm	11 x 38		6 x 45	10 x 60	12 x 75	12 x 82	12 x 100	14 x 100	17 x 100	24 x 100	34 x 100	44 x 100	
Cat. No.	2078	0536	Rhesus	-	0553 ²⁾	0501 ²⁾	0578 ²⁾	0500 ²⁾	0518 ²⁾	0519 ²⁾	0521 ²⁾	0526 ²⁾	cyto chambers ⁵⁾
 carrier Cat. No. 1425													
													
Cat. No.	1444		1432	1433	1438	1434		1431		1435	1436	1437	1452
boring Ø x L in mm	11.5 x 38		6.5 x 34	10.5 x 40	13.4 x 50	12.7 x 60		17.5 x 84		25.5 x 84	35.5 x 84	45.5 x 86	-
tubes per rotor	36		144	56	28	48		28		8	4	4	4
max. RCF ²⁾	3,885		3,969	3,913								3,801	2,739
radius in mm	139		142	140								136	98
run-up in sec	30												
run-down in sec, braked	32												
temperature in °C ¹⁾	-10												

capacity in ml	2.6–3.4	2.7–3	4–5.5	4.5–5	4.9	7.5–8.5	9–10	1.6–5	4–7	4–7	8.5–10	15	50	50	
Ø x L in mm	13 x 65	11 x 66	15 x 75	11 x 92	13 x 90	15 x 92	16 x 92	13 x 75	13 x 100	16 x 75	16 x 100	17 x 120	29 x 115	29 x 115	
Cat. No.	blood collection/urine tubes											0509	0513	-	
 carrier Cat. No. 1425															
															
Cat. No.	1438		1441	1438		1441	1439	1438		1441		1442	1443	1737	
boring Ø x L in mm	13.4 x 50		16.5 x 50	13.4 x 50		16.5 x 50	17 x 45	13.4 x 50		16.5 x 50		17x90	30 x 90	30 x 90	
tubes per rotor	28						16	28					4		
max. RCF ²⁾	3,913											4,081			
radius in mm	140											146			
run-up in sec	30														
run-down in sec, braked	32														
temperature in °C ¹⁾	-10														

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

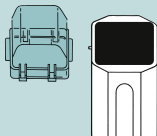
³⁾ When using these tubes, carrier 1427 cannot be closed with lid 1421.

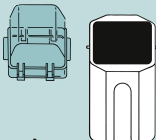
ROTORS AND ACCESSORIES



illustrated with carriers 1427

Swing-out rotor Cat. No. (without carriers) 1494

capacity in ml	1.5/2.0		1	3	5	6	7	9	15	25	50	1.1–1.4	2.6–3.4	2.7–3	4–5.5
Ø x L in mm	11 x 38		6 x 45	10 x 60	12 x 75	12 x 82	12 x 100	14 x 100	17 x 100	24 x 100	34 x 100	8 x 66	13 x 65	11 x 66	15 x 75
Cat. No.	2078/0536	-	Rhesus	-	0553 ²⁾	0501 ²⁾	0578 ²⁾	0500 ²⁾	0518 ²⁾	0519 ²⁾	0521 ²⁾	blood collection/urine tubes			
lid Cat. No. 1421  carrier Cat. No. 1427															
															
Cat. No.	5277	5273	1357	1327	1732	5229	5230	5237	5231	5232	5233	5278	1732	5229	5279
boring Ø x L in mm	11.5 x 38	11.2 x 38	6.5 x 23	10.5 x 23	13.4 x 58	12.5 x 44	12.4 x 87	15 x 73	17.8 x 87	25.5 x 87	35.5 x 87	9 x 41	13.4 x 58	12.5 x 44	15.6 x 41
tubes per rotor	36	48	120	48	32	48	48	24		8	4	48	32	48	20
max. RCF ²⁾	4,109	4,053	4,081	4,053	4,025	4,053	3,941					4,109	4,025	4,053	4,109
radius in mm	147	145	146	145	144	145	141					147	144	145	147
run-up in sec	30														
run-down in sec, braked	32														
temperature in °C ¹⁾	-7														

capacity in ml	4.5–5	4.9	7.5–8.5	9–10	1.6–5	4–7		8.5–10	15	50	25	30	50
Ø x L in mm	11x92	13x90	15x92	16x92	13x75	13x100	16x75	16x100	17x120	29x115	25x90	25x110	29x115
Cat. No.	blood collection/urine tubes								0509	0513	-	-	-
lid Cat. No. 1421  carrier Cat. No. 1427													
													
Cat. No.	5230	1732	5279	5271 ⁴⁾	1732	5271	5231	5275	5276	1731	5272	5233 ⁴⁾	
boring Ø x L in mm	12.4x87	13.4x58	15.6x41	17x66	13.4x58	17x66	17.8x87	17x90	30x90	25.5x72	25.5x85	35.5x87	
tubes per rotor	48	32	20		32	20	24	4		8	4	4	
max. RCF ²⁾	3,941	4,025	4,109	3,969	4,025	3,969	3,941	4,165	4,053	3,969	4,025	4,053	
radius in mm	141	144	147	142	144	142	141	149	145	142	144	145	
run-up in sec	30												
run-down in sec, braked	32												
temperature in °C ¹⁾	-7												

⁴⁾ Please remove the inserts.

⁵⁾ Please ask for our detailed cyto brochure.

ROTORS AND ACCESSORIES



illustrated with carriers 1495 and lids 1492
with bio-containment⁶⁾

Swing-out rotor Cat. No. (without carriers) **1494**

capacity in ml	1.5	2.0	1	3	5/6/7	9	15	25	50	85	100	1.1–1.4	2.7–5	2.6–4.9	4–8.5	
Ø x L in mm	11 x 38		6x45	10x60	12 x 75/82/100	14 x 100	17 x 100	24 x 100	34 x 100	38 x 101	44 x 100	8 x 66	11 x 66/92	13 x 65/90	15 x 75/92	
Cat. No.	2078	0536	Rhesus	-	0553 ²⁾ /0501 ²⁾ /0578 ²⁾	0500 ²⁾	0518 ²⁾	0519 ²⁾	0521 ²⁾	0549	0526 ²⁾	blood/urine tubes				
lid with bio-containment ⁶⁾ Cat. No. 1492																
																
carrier Cat. No. 1495																
Cat. No.	1351	1339	1343	1383	1329	1330	1331	1396	0761	1457	1383	1459				
boring Ø x L in mm	11.2 x 38	6.5 x 34	10.5 x 43	13.4 x 48	17.6 x 91	25.2 x 87	35.2 x 87	38.5 x 92	45.6 x 98	9 x 47	13.4 x 48	15.6 x 47				
tubes per rotor	20	108	36	20	16	4				28	20	16				
max. RCF ²⁾	3,857	4,081	4,137	4,025	3,997	3,829	3,801	4,109	4,025	4,053	4,025	4,053				
radius in mm	138	146	148	144	143	137	136	147	144	145	144	145				
run-up in sec	30															
run-down in sec, braked	32															
temperature in °C ¹⁾	-10															

capacity in ml	9 – 10	10	1.6 – 7	4 – 10	15	30	50	12	25	30	50	10	30	50	85	30	
Ø x L in mm	16 x 92	15 x 102	13 x 75/100	16 x 75/100	17 x 120	25 x 110	29 x 115	17 x 100	25 x 90	25 x 110	29 x 115	16 x 80	26 x 95	29 x 107	38 x 106	44 x 105	
Cat. No.	blood/urine tubes				0509	-	0513	-	-	-	-	-	0545	0546	0547	0534 ⁷⁾	
lid with bio-containment ⁶⁾ Cat. No. 1492																	chrome bath tube
																	
Cat. No.	1329	1329 ⁴⁾	1383	1348	1347	1365	1384	6311	1363	1365	6318	1348	4417	4416	1396	0765	
boring Ø x L in mm	17.6 x 91		13.4 x 48	16.5 x 56	17 x 90	26 x 80	30 x 90	17 x 80	26 x 72	26 x 80	29.5 x 80	16.5 x 56	26 x 83	29 x 93	38.5 x 92	45.9 x 98	
tubes per rotor	16		20	16	4							16	4			4	
max. RCF ²⁾	3,997	4,193	4,025		4,193				3,689	4,193		4,025	3,857	4,137	4,109	3,997	
radius in mm	143	150	144		150				132	150		144	138	148	147	143	
run-up in sec	30																
run-down in sec, braked	32																
temperature in °C ¹⁾	- 10																

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

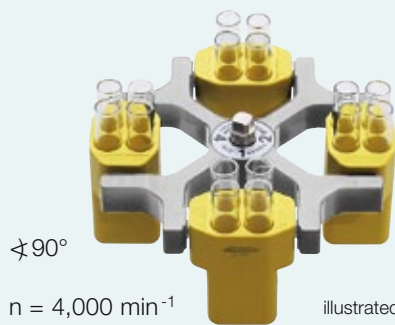
²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

⁴⁾ Please remove the inserts.

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

ROTORS AND ACCESSORIES

Swing-out rotor, 4-place

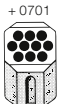


±90°

$n = 4,000 \text{ min}^{-1}$
max. RCF 2,665

illustrated with carriers 1369

Cat. No. (without carriers) **1624**

capacity in ml	5		6	7	9		15		20	25	45	50	1.1–1.4	2.6–3.4	4–5.5
Ø x L in mm	12 x 75		12 x 82	12 x 100	14 x 100		17 x 100		21 x 100	24 x 100	31 x 100	34 x 100	8 x 66	13 x 65	15 x 75
Cat. No.	0553 ²⁾		0501 ²⁾	0578 ²⁾	0500 ²⁾		0518 ²⁾		-	0519 ²⁾	-	0521 ²⁾	blood/urine tubes		
 rotor Cat. No. 1624															
	 with decanting aid		 with decanting aid			 + 0701							 + 0701	 + 0716	
Cat. No.	1369-91	1372	1369-92		1370	1741	1369	1742	1346	1745	1345	1746	1741	1742	
boring Ø x L in mm	11.5 x 64.4	13.5 x 65	11.5 x 71.5		14.6 x 74	14.6 x 78	17.6 x 74	17.6 x 78	21.5 x 74	26 x 78	32 x 74	35 x 78	14.6 x 78	17.6 x 78	
tubes per rotor	16	68	16		20	40	16	28	8		4		40	28	
max. RCF ²⁾	2,057	2,164	2,308			2,415	2,308	2,451	2,361	2,451	2,361	2,451	2,415	2,325	
radius in mm	115	121	129			135	129	137	132	137	132	137	135	130	
run-up in sec	20														
run-down in sec, braked	25														
temperature in °C ¹⁾	- 17				- 15	- 17	- 19	- 17	- 15	- 17	- 15		- 19		

capacity in ml	4.9	1.6–5	4–7	8.5–10	30	suspensions for cyto accessories	
Ø x L in mm	13 x 90	13 x 75	16 x 75	16 x 100	26 x 95		
Cat. No.	blood collection/urine tubes				0545	cyto chambers⁵⁾	
 rotor Cat. No. 1624							
Cat. No.	1741	1742	1369	1745	1660	1680	
boring Ø x L in mm	14.6 x 78	17.6 x 78	17.6 x 74	26 x 78	-		
tubes per rotor	40	28	16	8	4		
max. RCF ²⁾	2,451	2,325	2,308	2,451	1,646	1,467	
radius in mm	137	130	129	137	92	82	
run-up in sec	20						
run-down in sec, braked	25						
temperature in °C ¹⁾	-15	-19	-17	-15	-16	-17	

capacity in ml	1.5	2.0	1	3	4
Ø x L in mm	11 x 38		6 x 45	10 x 60	12 x 60
Cat. No.	2078	0536	Rhesus	-	-
 carrier Cat. No. 1366					
					
Cat. No.	5277		1357	1327	1326
boring Ø x L in mm	11.5 x 38		6.5 x 23	10.5 x 23	12.5 x 44
tubes per rotor	36		120	48	
max. RCF ²⁾	1,968		1,950	1,932	
radius in mm	110		109	108	
run-up in sec	20				
run-down in sec, braked	25				
temperature in °C ¹⁾	- 17				

⁵⁾ Please ask for our detailed cyto brochure.

⁷⁾ A rubber stopper for closing the tube is available under Cat. No. 0535.

The tube may not be centrifuged with the stopper. The stopper only serves to close the tube for agitating or mixing.

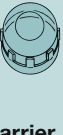
ROTORS AND ACCESSORIES



illustrated with carriers 1481 and lids 1492
with bio-containment⁶⁾

Swing-out rotor Cat. No. (without carriers) 1624

capacity in ml	1.5	2.0	1	3	5/6/7	9	15	25	50	85	100	1.1–1.4	2.7–5	2.6–4.9	4–8.5
Ø x L in mm	11 x 38	6 x 45	10 x 60	12 x 75/82/100	14 x 100	17 x 100	24 x 100	34 x 100	38 x 101	44 x 100	8 x 66	11 x 66/92	13 x 65/90	15 x 75/92	
Cat. No.	2078	0536	Rhesus	-	0553 ²⁾ /0501 ²⁾ /0578 ²⁾	0500 ²⁾	0518 ²⁾	0519 ²⁾	0521 ²⁾	0549	0526 ²⁾	blood/urine tubes			
lid with bio-containment ⁶⁾ Cat. No. 1492  carrier Cat. No. 1481															
															
Cat. No.	1351	1339	1343	1383	1329	1330	1331	1396	0761	1457	1383	1459			
boring Ø x L in mm	11.2 x 38	6.5x34	10.5 x 43	13.4 x 48	17.6x91	25.2x87	35.2x87	38.5 x 92	45.6 x 98	9 x 47	13.4 x 48	15.6 x 47			
tubes per rotor	20	108	36	20	16	4				28	20	16			
max. RCF ²⁾	2,451	2,594	2,630	2,558	2,540	2,433	2,415	2,612	2,558	2,576	2,558	2,576			
radius in mm	137	145	147	143	142	136	135	146	143	144	143	144			
run-up in sec	20														
run-down in sec, braked	25														
temperature in °C ¹⁾	-15														

capacity in ml	9–10	10	1.6–7	4–10	15	50	12	25	30	50	10	30	50	85	30
Ø x L in mm	16 x 92	15 x 102	13 x 75/100	16 x 75/100	17 x 120	29 x 115	17 x 100	25 x 90	25 x 110	29 x 115	16 x 80	26 x 95	29 x 107	38 x 106	44 x 105
Cat. No.	blood/urine tubes				0509	0513	-	-	-	-	-	0545	0546	0547	0534 ⁷⁾
lid with bio-containment ⁶⁾ Cat. No. 1492  carrier Cat. No. 1481															
															
Cat. No.	1329	1329 ⁴⁾	1383	1348	1347	1384	6311	1363	1365	6318	1348	4417	4416	1396	0765
boring Ø x L in mm	17.6 x 91	13.4 x 48	16.5 x 56	17 x 90	30 x 90	17 x 80	26 x 72	26 x 80	29.5 x 80	16.5 x 56	26 x 83	29 x 93	38.5 x 92	45.9 x 98	
tubes per rotor	16	20	16	4				16			4				
max. RCF ²⁾	2,540	2,665	2,558	2,522	2,665	2,343	2,665	2,522	2,451	2,630	2,612	2,540			
radius in mm	142	149	143	141	149	131	149	141	137	147	146	142			
run-up in sec	20														
run-down in sec, braked	25														
temperature in °C ¹⁾	-15														

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

⁴⁾ Please remove the inserts.

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

ROTORS AND ACCESSORIES

Swing-out rotor, 4-place

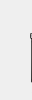

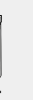


$n = 4,500 \text{ min}^{-1}$
max. RCF 3,328

illustrated with carriers 1490 and lids 1492
with bio-containment⁶⁾

Cat. No. (without carriers) **1324**

capacity in ml	1.5	2.0	1	3	4	5	6	7	9	15	25	50	85	100
Ø x L in mm	11 x 38		6 x 45	10 x 60	10 x 88	12 x 75	12 x 82	12 x 100	14 x 100	17 x 100	24 x 100	34 x 100	38 x 101	44 x 100
Cat. No.	2078	0536	Rhesus	-	-	0553 ²⁾	0501 ²⁾	0578 ²⁾	0500 ²⁾	0518 ²⁾	0519 ²⁾	0521 ²⁾	0549	0526 ²⁾
lid with bio-containment ⁶⁾ Cat. No. 1492   carrier Cat. No. 1490														
														
Cat. No.	1351		1339	1343		1383			1329		1330	1331	1396	0761
boring Ø x L in mm	11.2 x 38		6.5 x 34	10.5 x 43		13.4 x 48			17.6 x 91		25.2 x 87	35.2 x 87	38.5 x 92	45.9 x 100.5
tubes per rotor	20		108	36		20			16		4			
max. RCF ²⁾	3,056		3,237	3,283		3,192			3,170		3,034	3,011	3,260	3,192
radius in mm	135		143	145		141			140		134	133	144	141
run-up in sec	27													
run-down in sec, braked	30													
temperature in °C ¹⁾	-6													

capacity in ml	1.1–1.4	2.6–3.4	4.9	2.7–3	4.5–5	4–5.5	7.5–8.5	9–10	10	1.6–5	4–7	4–7	8.5–10
Ø x L in mm	8x66	13x65	13x90	11 x 66	11 x 92	15 x 75	15 x 92	16 x 92	15 x 102	13 x 75	13 x 100	16x75	16x100
Cat. No.	blood collection/urine tubes												
lid with bio-containment ⁶⁾ Cat. No. 1492   carrier Cat. No. 1490													
													
Cat. No.	1457	1383				1459		1329 1329 ⁴⁾		1383		1348	
boring Ø x L in mm	9x47	13.4x48				15.6x47		17.6x91		13.4x48		16.5x56	
tubes per rotor	28	20				16				20		16	
max. RCF ²⁾	3,215	3,192				3,215		3,170	3,328	3,192		3,147	
radius in mm	142	141				142		140	147	141		139	
run-up in sec	27												
run-down in sec, braked	30												
temperature in °C ¹⁾	-6												

⁷⁾ A rubber stopper for closing the tube is available under Cat. No. 0535.
The tube may not be centrifuged with the stopper. The stopper only serves to close the tube for agitating or mixing.

ROTORS AND ACCESSORIES



illustrated with carriers 1398

Swing-out rotor 1324

capacity in ml	15		30	50	12	25	30	50	10	30	50	85	30
Ø x L in mm	17 x 120		25 x 110	29 x 115	17 x 100	25 x 90	25 x 110	29 x 115	16 x 80	26 x 95	29 x 107	38 x 106	44 x 105
Cat. No.	0509		-	0513	-	-	-	-	-	0545	0546	0547	0534 ⁷⁾
lid with bio-containment ⁶⁾ Cat. No. 1492  carrier Cat. No. 1490 		 ³⁾											 chrome bath tube
													 rubber insert
Cat. No.	1347	1356	1365	1384	6311	1363	1365	6318	1348	4417	4416	1396	0765
boring Ø x L in mm	17 x 90	17 x 107	26 x 80	30 x 90	17 x 80	26 x 72	26 x 80	29.5 x 80	16.5 x 56	26 x 83	29 x 93	38.5 x 92	45.9 x 100.5
tubes per rotor	4	12	4						16	4			
max. RCF ²⁾	3,328		3,328			2,920	3,328		3,147	3,056	3,283	3,260	3,170
radius in mm	147		147			129	147		139	135	145	144	140
run-up in sec	27												
run-down in sec, braked	30												
temperature in °C ¹⁾	- 6												

capacity in ml	5	6	7	9	15	2.6–3.4	2.7–3	4–5.5	4.5–5	4.9	9–10	10	1.6–5	4–7	4–7	12	
Ø x L in mm	12 x 75	12 x 82	12 x 100	14 x 100	17 x 100	13 x 65	11 x 66	15 x 75	11 x 92	13 x 90	16 x 92	15 x 102	13 x 75	13 x 100	16 x 75	17 x 102	
Cat. No.	0553 ²⁾	0501 ²⁾	0578 ²⁾	0500 ²⁾	0518 ²⁾	blood collection/urine tubes											
 carrier Cat. No. 1398																	
																	
Cat. No.	1486			1482A		1486		1482A		1486		1482A		1486		1487	
boring Ø x L in mm	13.4 x 57.5			17.5 x 81		13.4 x 57.5		17.5 x 81		13.4 x 57.5		17.5 x 81		13.4 x 57.5		17.5 x 74	
tubes per rotor	20			16		20		16		20		16		20		12	
max. RCF ²⁾	3,147			3,192		3,147		2,875		3,147		3,192		3,147		3,102	
radius in mm	139			141		139		127		139		141		139		137	
run-up in sec	27																
run-down in sec, braked	30																
temperature in °C ¹⁾	- 6																

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

³⁾ When using these tubes, carrier 1490 cannot be closed with lid 1492.

ROTORS AND ACCESSORIES

Swing-out rotor, 6-place



✂ 90°

n = 4,000 min⁻¹
max. RCF 2,701

Cat. No. (with carriers) 1619

Swing-out rotor, 8-place



✂ 45°

n = 5,000 min⁻¹
max. RCF 3,857

Cat. No. (with carriers) 1617

Swing-out rotor, 12-place



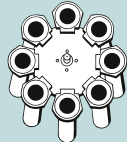
✂ 55°
(with carriers 1127)
✂ 60°
(with carriers 1122)
✂ 80°
(with carriers 1621)

n = 5,000 min⁻¹
max. RCF 4,193

Cat. No. (without carriers) 1628

capacity in ml	15	50	12	50
Ø x L in mm	17 x 120	29 x 115	17 x 100	29 x 115
Cat. No.	0509	0513	-	-
 carrier Cat. No. 1398				
				
Cat. No.	1483A	1484	1482A	1484⁴⁾
boring Ø x L in mm	17 x 100	30 x 98	17.5 x 81	30 x 98
tubes per rotor	16	4	16	4
max. RCF ²⁾	3,305	3,260	3,192	3,260
radius in mm	146	144	141	144
run-up in sec	27			
run-down in sec, braked	30			
temperature in °C ¹⁾	-6			

capacity in ml	15	50
Ø x L in mm	17 x 120	29 x 115
Cat. No.	0509	0513
 rotor Cat. No. 1619		
		
Cat. No.	1462-A	-
boring Ø x L in mm	17 x 84	30 x 87.5
tubes per rotor	6	
max. RCF ²⁾	2,701	
radius in mm	151	
run-up in sec	20	
run-down in sec, braked	22	
temperature in °C ¹⁾	-15	

capacity in ml	15	50
Ø x L in mm	17 x 120	29 x 115
Cat. No.	0509	0513
 rotor Cat. No. 1617		
		
Cat. No.	1462-A	-
boring Ø x L in mm	17 x 84	30 x 94.5
tubes per rotor	8	
max. RCF ²⁾	3,857	
radius in mm	138	
run-up in sec	20	
run-down in sec, braked	19	
temperature in °C ¹⁾	-10	

capacity in ml	5	10	15	2.6–3.4	2.7–3	4–5.5	7.5–8.5	1.6–5	4–7	8.5–10
Ø x L in mm	12 x 75	17x70	17 x 100	13x65	11 x 66	15x75	15x92	13x75	16x75	16 x 100
Cat. No.	0553 ²⁾	2079	0518 ²⁾	blood collection/urine tubes						
 rotor Cat. No. 1628										
										
Cat. No.	1127	1122	1621	1127	1122	1621	1127	1122	1621	
boring Ø x L in mm	13.2 x 53	17.5 x 53	17.5 x 79	13.2 x 53	17.5 x 53	17.5 x 79	13.2 x 53	17.5 x 53	17.5 x 79	
tubes per rotor	12									
max. RCF ²⁾	3,494	3,522	4,193	3,494	3,522	4,193	3,494	3,522	4,193	
radius in mm	125	126	150	125	126	150	125	126	150	
run-up in sec	16									
run-down in sec, braked	16									
temperature in °C ¹⁾	- 15	- 10	- 15	- 15	- 10	- 15	- 15	- 15	- 10	

⁴⁾ Please remove the inserts.

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

⁷⁾ A rubber stopper for closing the tube is available under Cat. No. 0535.

The tube may not be centrifuged with the stopper. The stopper only serves to close the tube for agitating or mixing.

ROTORS AND ACCESSORIES

Swing-out rotor, 8-place

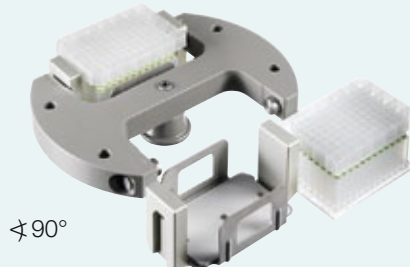


±90°

$n = 4,000 \text{ min}^{-1}$
max. RCF 2,415

Cat. No. (without carriers) **1611**

Swing-out rotor, 2-place

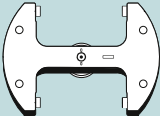
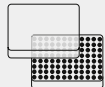


±90°

$n = 4,000 \text{ min}^{-1}$
max. RCF 2,218

Cat. No. (without carriers) **1460**

capacity in ml	5	6	7	10		15	2.6–3.4	2.7–3	4–5.5	4.5–5	7.5–8.5	1.6–5	4–7	4–7	8.5–10
Ø x L in mm	12 x 75	12 x 82	12 x 100	17 x 70	13 x 100	17 x 100	13 x 65	11 x 66	15 x 75	11 x 92	15 x 92	13 x 75	13 x 100	16 x 75	16 x 100
Cat. No.	0553 ²⁾	0501 ²⁾	0578 ²⁾	2079	-	0518 ²⁾	blood collection/urine tubes								
 rotor Cat. No. 1611															
															
Cat. No.	1131		1643	1132	1643	1644	1131		1132	1643	1644	1131	1643	1132	1644
boring Ø x L in mm	13 x 53		13.2 x 81	17.5 x 53	13.2 x 81	17.5 x 81	13 x 53		17.5 x 53	13.2 x 81	17.5 x 81	13 x 53	13.2 x 81	17.5 x 53	17.5 x 81
tubes per rotor	8														
max. RCF ²⁾	1,914		2,415	1,914	2,415		1,914		2,415		1,914	2,415	1,914	2,415	
radius in mm	107		135	107	135		107		135		107	135	107	135	
run-up in sec	20														
run-down in sec, braked	20														
temperature in °C ¹⁾	- 16														

D x W x H in mm	86 x 128 x 15/17.5		86 x 128 x 22	86 x 128 x 44.5	86 x 128 x 46	86 x 128 x 83	59 x 84 x 11	82 x 124 x 20	-
capacity in ml	-								0.2
Cat. No.	MTP		CP	DWP	MS	QP	Microtest plate (Terasaki)	PCR plate, 96 wells	PCR strips
 rotor Cat. No. 1460									
								 	
Cat. No.	1453-A							1453-A + 1485	
boring Ø x L in mm	-							-	
plates/strips pro rotor	10	8	6	2	2	2	4	2	24 x 8
max. RCF ²⁾	2,218								
radius in mm	124								
run-up in sec	40								
run-down in sec, braked	45								
temperature in °C ¹⁾	- 6								

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

⁸⁾ For centrifugation at high speeds, we recommend to use form-fitting, phenol-resistant adapters 2031.

ROTORS AND ACCESSORIES

Swing-out rotor, 4-place



±90°

n = 4,000 min⁻¹
max. RCF 1,896

Cat. No. (without carriers) **1645**

Angle rotor, 12-place



±45°

n = 15,000 min⁻¹
max. RCF 16,602

Cat. No. **1612**

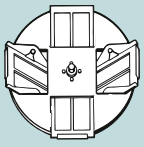
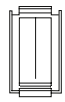
Angle rotor, 24-place



±40°

n = 15,000 min⁻¹
max. RCF 21,382

Cat. No. **1420-B**

D x W x H in mm	86 x 128 x 15/17.5		86 x 128 x 22
Cat. No.	MTP		CP
			
rotor Cat. No. 1645			
Cat. No.	4345		
boring Ø x L in mm	-		
plates per rotor	4		
max. RCF ²⁾	1,896		
radius in mm	106		
run-up in sec	20		
run-down in sec, braked	19		
temperature in °C ¹⁾	- 16		

capacity in ml	0.2	0.4	0.5	0.8	1.5	2.0
Ø x L in mm	6 x 18	6 x 45	8 x 30	8 x 45	11 x 38	
Cat. No.	-	-	-	-	2078	0536
 rotor Cat. No. 1612 lid incl.						
						
Cat. No.	2024		2023		2031 ⁶⁾	-
boring Ø x L in mm	6 x 40		8 x 40		10.2x19	11.4 x 39
tubes per rotor	12					
max. RCF ²⁾	16,602					
radius in mm	66					
run-up in sec	25					
run-down in sec, braked	23					
temperature in °C ¹⁾	- 10					

capacity in ml	0.2	0.4	0.5	0.8	1.5	2.0
Ø x L in mm	6 x 18	6 x 45	8 x 30	8 x 45	11 x 38	
Cat. No.	-	-	-	-	2078	0536
 rotor Cat. No. 1420-B lid incl.						
						
Cat. No.	2024		2023		2031 ⁸⁾	-
boring Ø x L in mm	6 x 40		8 x 40		10.2 x 19	11.2 x 42.6
tubes per rotor	24					
max. RCF ²⁾	21,382					
radius in mm	85					
run-up in sec	25					
run-down in sec, braked	24					
temperature in °C ¹⁾	- 4					

Optional for rotor 1420-B:



lid **Cat. No. 2425**
with bio-containment⁶⁾,
autoclavable



lid **Cat. No. 2423**
with bio-containment⁶⁾,
autoclavable, phenol-resistant

ROTORS AND ACCESSORIES

Angle rotor, 30-place



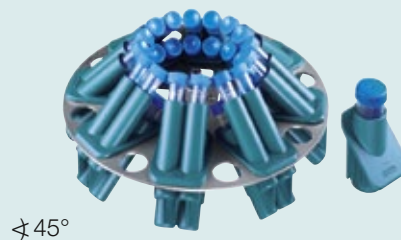
45°

n = 14,000 min⁻¹
max. RCF 21,255

with bio-containment⁶⁾

Cat. No. 1689-A

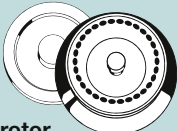
Angle rotor, 8-place



45°

n = 4,500 min⁻¹
max. RCF 3,283

Cat. No. (without adapters) 1418

capacity in ml	0.2	0.4	0.5	0.8	1.5	2.0
Ø x L in mm	6 x 18	6 x 45	8 x 30	8 x 45	11 x 38	
Cat. No.	-	-	-	-	2078	0536
lid with bio-containment ⁶⁾ incl.  rotor Cat. No. 1689-A						
						
Cat. No.	2024		2023		2031 ⁶⁾	-
boring Ø x L in mm	6 x 40		8 x 40		10.2 x 19	11.2 x 40.9
tubes per rotor	30					
max. RCF ²⁾	21,255					
radius in mm	97					
run-up in sec	35					
run-down in sec, braked	31					
temperature in °C ¹⁾	0					

capacity in ml	5	15	1.1–1.4	2.6–3.4	2.7–3	9–10	1.6–5	12	15	50	12	50	50
Ø x L in mm	12 x 75	17 x 100	8 x 66	13 x 65	11 x 66	16 x 92	13 x 75	17 x 102	17 x 120	29 x 115	17 x 100	29 x 115	29 x 107
Cat. No.	0553 ²⁾	0518 ²⁾	blood collection/urine tubes						0509	0513	-	-	0546
 rotor Cat. No. 1418													
	+ 1054-A 	+ 0716 	+ 1054-A 			+ 0716 	+ 1054-A 	+ 0716 	+ E 2109 	+ E 2110-A 	+ 0716 		
Cat. No.	1467	1467	1467			1467	1467	1467	1467	1468	1467	1468	
boring Ø x L in mm	17.4 x 91	17.4 x 91	17.4 x 91			17.4 x 91	17.4 x 91	17.4 x 91	17.4 x 91	30.2 x 92	17.4 x 91	30.2 x 92	
tubes per rotor	32	32	32			32	32	32	32	8	32	8	
max. RCF ²⁾	2,762	3,215	2,762			3,215	2,762	3,215	3,283	3,147	3,215	3,147	
radius in mm	122	142	122			142	122	142	145	139	142	139	
run-up in sec	30												
run-down in sec, braked	31												
temperature in °C ¹⁾	- 11												

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000.

⁶⁾ Tested by the TÜV in conformity with DIN EN 61010, section 2-020.

⁹⁾ 6-place adapter set for conical 50 ml tubes with screw cap.

ROTORS AND ACCESSORIES

Angle rotor, 6-place



✱35°

n = 9,000 min⁻¹
max. RCF 9,509

Cat. No. 1620A

Angle rotor, 12-place



✱35°

n = 6,000 min⁻¹
max. RCF 4,146

Cat. No. 1613

capacity in ml	1.5	2.0	15	50	75	85	7.5–8.5	9–10	10	8.5–10	15	50	10	30	50	85	
Ø x L in mm	11 x 38		17 x 100	34 x 100	35 x 105	38 x 101	15 x 92	16 x 92	15 x 102	16 x 100	17 x 120	29 x 115	16 x 80	26 x 95	29 x 107	38 x 106	
Cat. No.	2078	0536	0518 ²⁾	0521 ²⁾	0548	0549	blood/urine tubes				0509	0513	-	0545	0546	0547	
 rotor Cat. No. 1620A																	
																	
Cat. No.	1449		1451	1463		-	1451				1466	1454	1646 ⁹⁾	1448	1447	1446	-
boring Ø x L in mm	11.4 x 39		17.5 x 92	35 x 89		38.6 x 90.2	17.5 x 92				17 x 106	29.8 x 97	-	16.5 x 74	26 x 85	29 x 92	38.6 x 90.2
tubes per rotor	24		6	6		6	6				6	3	6	12	6		
max. RCF ²⁾	9,237		8,784	9,327		9,509	8,784				8,965			8,784	8,603	9,056	9,509
radius in mm	102		97	103		105	97				99			97	95	100	105
run-up in sec	30																
run-down in sec, braked	30																
temperature in °C ¹⁾	- 11																

capacity in ml	4	5	6	15	1.1–1.4	2.6–3.4	2.7–3	4.5–5	4.9	7.5 – 10	10	1.6 – 5	4 – 7	8 – 10	15		
Ø x L in mm	10x88	12x75	12x82	17x100	8x66	13x65	11x66	11x92	13x90	15/16x92	15x102	13x75	13x100	16x100/125	17x120		
Cat. No.	-	0553 ²⁾	0501 ²⁾	0518 ²⁾	blood collection/urine tubes										0509		
 rotor Cat. No. 1613																	
																	
Cat. No.	6305	1054-A	-	1054-A	-							1054-A	1058	-			
boring Ø x L in mm	11.5x67.5	13.5 x 60	17.7 x 88	13.5x60	17.7 x 88							13.5x60	13.5 x 79	17.7 x 88			
tubes per rotor	12												12	12	12	6	6
max. RCF ²⁾	3,502	3,300	4,146	3,300	4,146							3,300	4,146				
radius in mm	87	82	103	82	103							82	103				
run-up in sec	15																
run-down in sec, braked	15																
temperature in °C ¹⁾	- 16																

ROTORS AND ACCESSORIES

Angle rotor, 12-place



✂ 35°

n = 12,000 min⁻¹
max. RCF 16,582

Cat. No. 1615

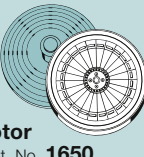
Hematocrit rotor, 24-place



n = 15,000 min⁻¹
max. RCF 21,382

Cat. No. 1650

capacity in ml	4	5	6	15	1.1–1.4	2.6–3.4	2.7–3	4.5–5	4.9	7.5–10	10	1.6–5	4–7	8–10	15	
Ø x L in mm	10x88	12x75	12x82	17 x 100	8x66	13x65	11x66	11 x 92	13x90	15/16x92	15x102	13x75	13x100	16x100/125	17x120	
Cat. No.	-	0553 ²⁾	0501 ²⁾	0518 ²⁾	blood collection/urine tubes										0509	
 rotor Cat. No. 1615																
																
Cat. No.	6305	1054-A	-	1054-A	-							1054-A	1058	-	1647 ¹⁰⁾	
boring Ø x L in mm	11.5x67.5	13.5 x 60	17.7x88	13.5x60	17.7x88							13.5x60	13.5x79	17.7x88	17x104	
tubes per rotor	12											12	12	12	6	6
max. RCF ²⁾	14,006	13,201	16,582	13,201	16,582							13,201	16,582	16,582	15,455	
radius in mm	87	82	103	82	103							82	103	103	96	
run-up in sec	40															
run-down in sec, braked	40															
temperature in °C ¹⁾	- 2															

standard capillaries, heparinised	basic	self-sealing and mylar-coated
Cat. No.	2074	1071
 rotor Cat. No. 1650	 	
Cat. No.	2077	-
boring Ø x L in mm	-	
capillaries per rotor	24	
max. RCF ²⁾	21,382	
radius in mm	85	
run-up in sec	13	
run-down in sec, braked	9	
temperature in °C ¹⁾	- 1	

¹⁾ Lowest attainable temperature in precooled refrigerated centrifuges at max. speed. Lower temperatures can be reached by reducing the speed.

²⁾ Please note that the RCF values indicated refer only to rotor performance. The max. permissible RCF of tubes used should be verified with the individual manufacturers. The max. RCF for glass tubes annotated with footnote ²⁾ is 4,000

¹⁰⁾ Set of 6 carriers. For accommodating conical 15 ml tubes with screw cap, remove 6 opposite carriers from rotor No. 1615 and replace them with adapters No. 1647.

ROTORS FOR CYTOLOGY

Rotors **No. 1626** and **No. 1648** were especially designed for use in cytology. They are recommended for laboratories working exclusively in this area.



Cyto rotor, 6-place

$n = 4,000 \text{ min}^{-1}$
max. RCF 2,039

illustrated with suspensions 1660
and lids 1661 for cyto chambers

Cat. No. 1626 (without suspensions)



Cyto rotor, 8-place

$n = 4,000 \text{ min}^{-1}$
max. RCF 2,218

illustrated with suspensions 1680
for angle chambers

Cat. No. 1648 (with suspensions)

THE HETTICH CYTO SYSTEM

The Hettich Cyto System provides cytologists with an extensive range of accessories to prepare high-quality cyto sediments from the most diverse body fluids.

Benefits of the Hettich Cyto System include:

- Low purchase costs
- Easy assembly
- Valuable time savings
- Optimal cell yield
- Targeted preparation of dry or wet sediments
- First-class customer support

Detailed information on our comprehensive range of cyto accessories is provided in a separate brochure. We will be pleased to send one to you on request.

Application notes from our R&D department on the pre-paration of slides for the cytological examination of various body fluids can be obtained free of charge or downloaded from our homepage www.hettichlab.com

The assembly of a complete cyto insert is really easy



suspension 1660
for cyto chambers

suspension 1680
for angle chambers

TECHNOLOGY	UNIVERSAL 320		UNIVERSAL 320 R		
Benchtop centrifuge, without rotor	classic		cooled		
Power supply ^{*)}	200 – 240 V 1 ~	100 – 127 V 1 ~	200-240 V 1 ~	240 V 1 ~	115 – 127 V 1 ~
Frequency	50 – 60 Hz	50 – 60 Hz	50 Hz	60 Hz	60 Hz
Consumption	400 VA		800 VA		950 VA
Emission, Immunity	EN / IEC 61326 - 1, class B	FCC class B	EN / IEC 61326 - 1, class B		FCC class B
 Max. capacity	4 x 100 ml / 32 x 15 ml				
 Max. RPM (speed)	15,000 min ⁻¹				
Max. RCF	21,382				
Running time	1 sec – 99 min : 59 sec, ∞ continuous run, short cycle mode (impulse key)				
Dimensions (H x W x D)	346 x 395 x 520 mm		346 x 401 x 695 mm		
Weight	approx. 31 kg		approx. 53 kg		
Refrigeration					
Temperature control, infinitely variable	-		from - 20 to + 40 °C		
Cat. No.	1401	1401-01	1406	1406-01	

^{*)} Other voltages on request.



Hettich centrifuges comply with all relevant EU standards in effect and conform to the European level of quality and safety for medical devices. Evidence is provided by national and international test marks such as IEC 61010 or the CE conformity. The ISO 9001, ISO 13485 and ISO 14001 certificates accredited to the company bear witness to the extreme care and responsibility Hettich puts into the manufacturing of centrifuges and their accessories.

Andreas Hettich GmbH & Co. KG

Föhrenstr. 12
D - 78532 Tuttlingen
Germany
www.hettichlab.com
info@hettichlab.com
service@hettichlab.com

Phone +49 (0)7461 / 705 -0
Fax +49 (0)7461 / 705 -122

National Sales: -200
International Sales: -201
National Service: -202
International Service: -203



LAB TECHNOLOGY